Document No.	05SDIPD00311	Secret Level	<input type="checkbox"/> General Lv1 <input checked="" type="checkbox"/> Secrecy Lv2 <input type="checkbox"/> Confidential Lv3	
Version	A-2	Doc. Type	DHF	
Code	RP006-OD2	Code Model	ASUS-IPD-V100	
Scope	Endoscopic Video Image Processor (FDA)			
<h2>Instruction for Use of Endoscopic Video Image Processor</h2>				
Drafted by		Date		
Reviewed by		Date		
Reviewed by		Date		
Approved by		Date		
Date of Publication		Seal		
D i s t r i b u t i o n	<input type="checkbox"/> Quality Assurance Department	<input type="checkbox"/> Quality Inspection Department		
	<input type="checkbox"/> Laboratory	<input type="checkbox"/> Production Department		
	<input type="checkbox"/> R&D Department	<input type="checkbox"/> NPI Department		
	<input type="checkbox"/> Engineering Department	<input type="checkbox"/> PE Technical Department		
	<input type="checkbox"/> Planning Department	<input type="checkbox"/> Procurement Department		
	<input type="checkbox"/> Regulatory Affairs Department	<input type="checkbox"/> After-Sales Department		
	<input type="checkbox"/> Domestic Sales Department	<input type="checkbox"/> Domestic Marketing Department		
	<input type="checkbox"/> Overseas Sales Department	<input type="checkbox"/> Overseas Marketing Department		
	<input type="checkbox"/> Clinical Medical Product Department	<input type="checkbox"/> Administration & HR Department		

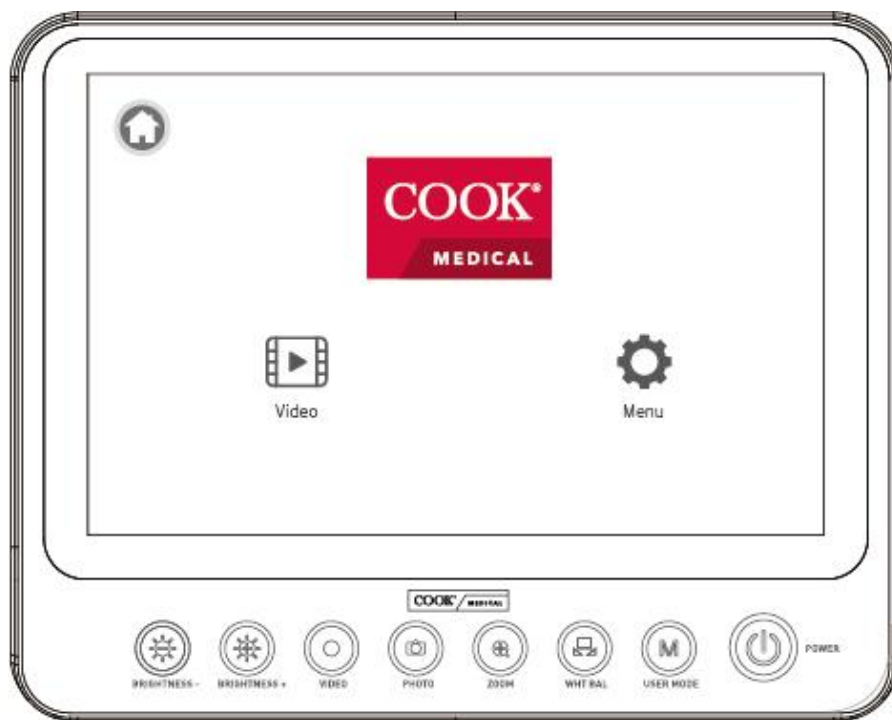
Change History

Rev.	Date	Revised contents	Prepared by	Change No.
A-0	2023.12.27	Initial release.	Xueping Zeng	20231227001
A-1	2024.03.12	Modify distributor addresses and optimize individual text descriptions.	Xueping Zeng	2020240312001
A-2	2024.10.18	Add information on Australian sponsor according to the suggestion of TGA audit	Yanting Liu	






User Manual

Endoscopic Video Image Processor

Model: ASUS-IPD-V100

Cook Medical LLC
400 Daniels Way
Bloomington, IN 47404, USA
Tel No. 800.457.4500
www.cookmedical.com

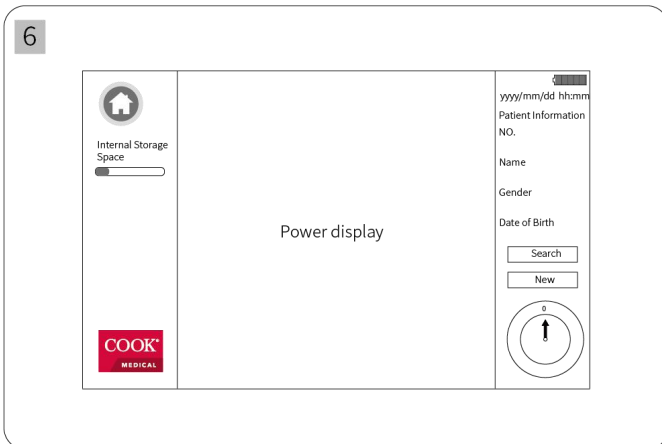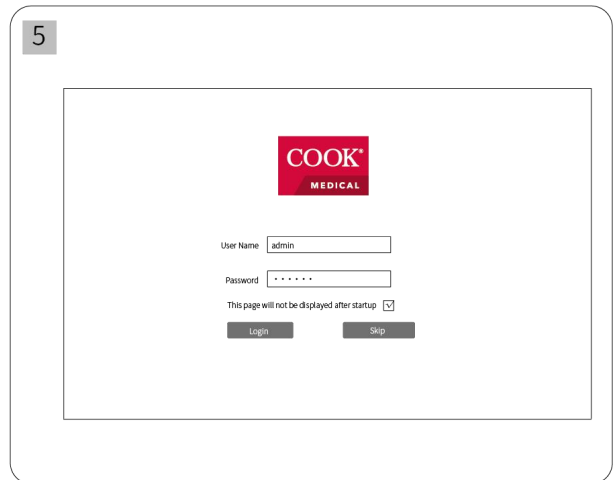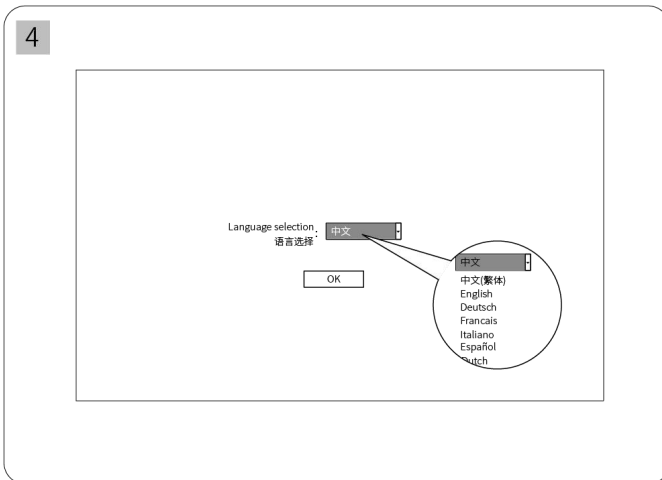
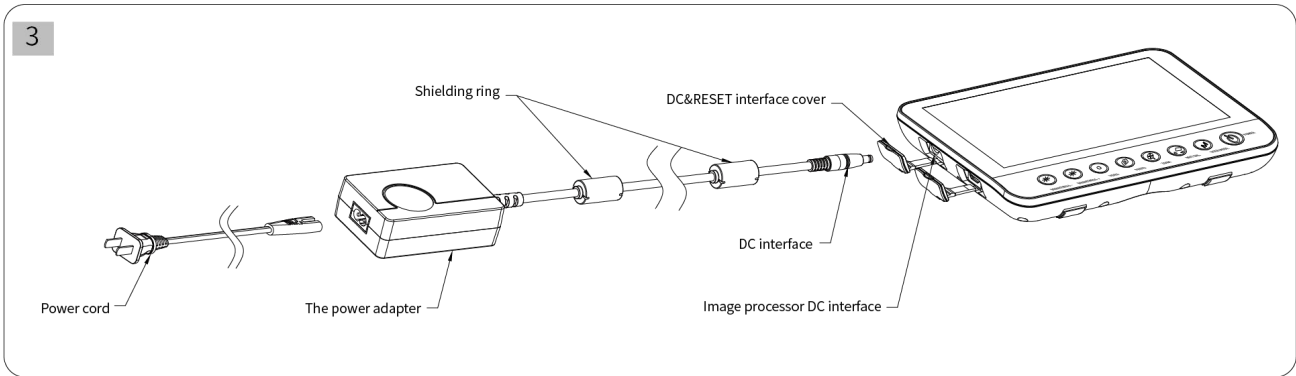
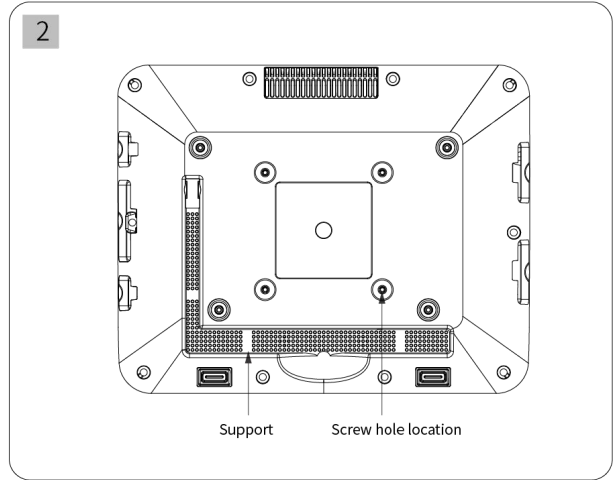
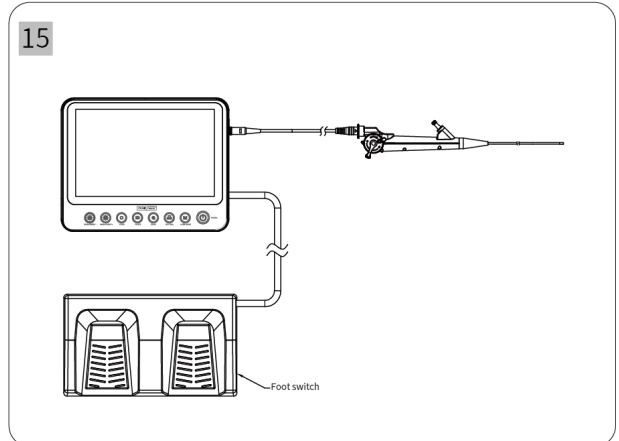
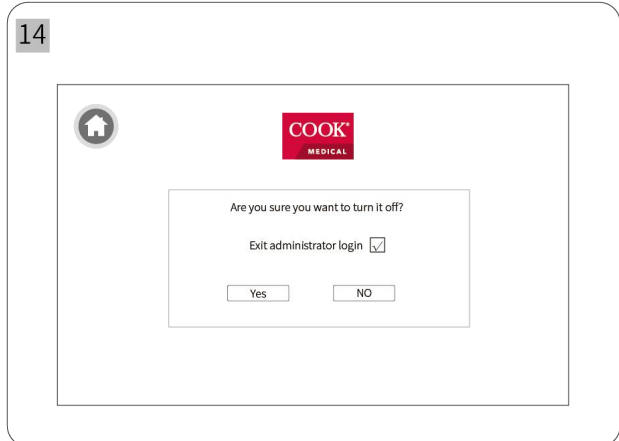
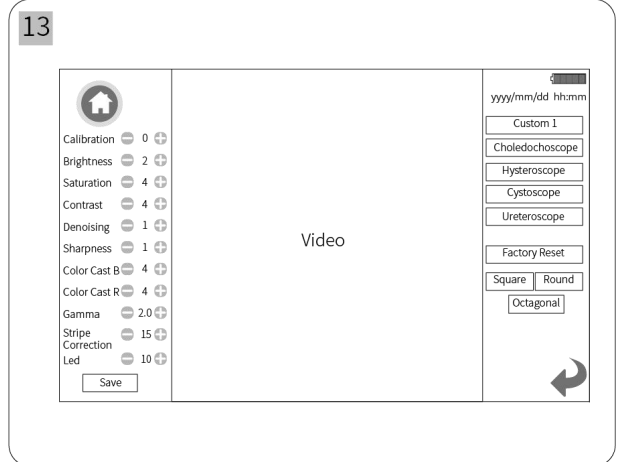
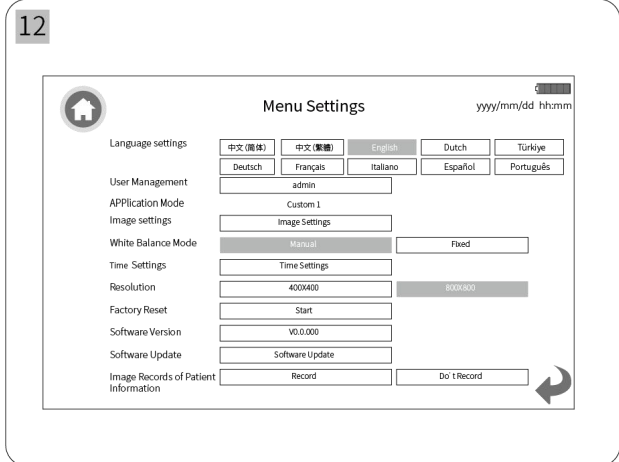
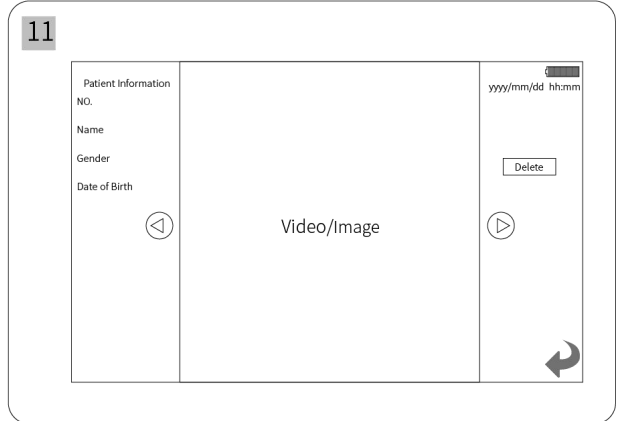
Table of Contents

1 User Guide	4
2 Safety Guidelines	5
3 Device Description	8
4 Control Operation of the Device Panel	8
5 Basic Operation	9
6 List of Supplies	12
7 Precautions for Use	13
8 Cleaning, Disinfection and Maintenance	14
9 Troubleshooting	14
10 Guarantee	15
11 Discard	15
12 Basic Parameters	16
13 Software-related Instructions	17
14 Attached Table	17

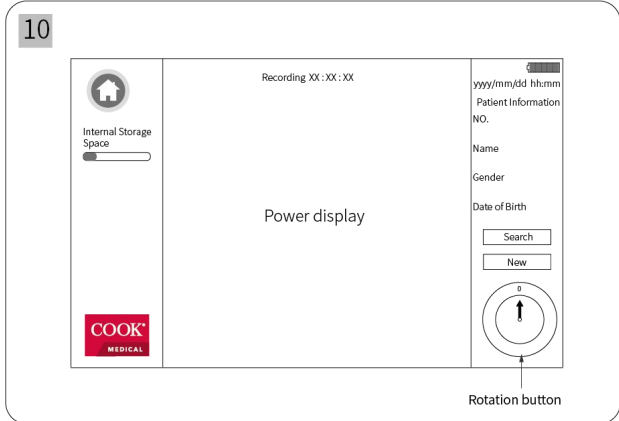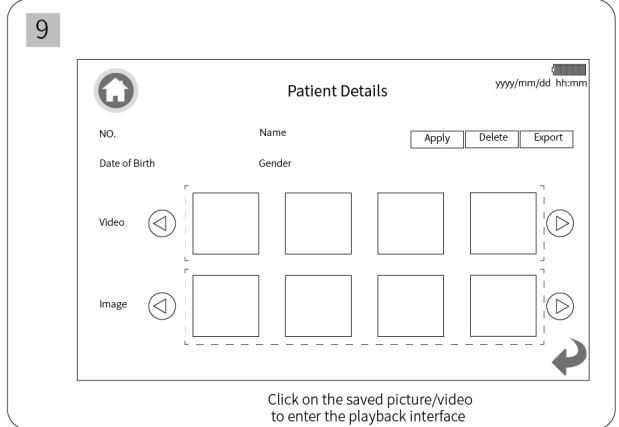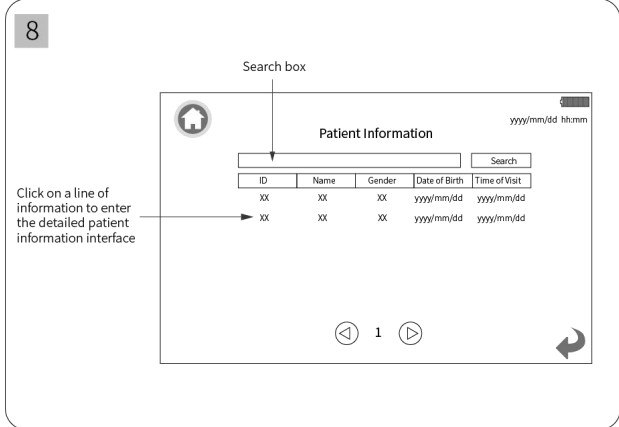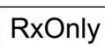
1 User Guide

1.1 Basic guide

- (1) This User Manual contains important contents that must be understood in the safe and correct use of this device;
- (2) This User Manual is part of the device. Therefore, during the entire service life of this device, the User Manual must be kept at the operating location of the device at any time;
- (3) This device must be operated by trained personnel with relevant experience. All personnel must read this User Manual carefully before using this device;
- (4) This User Manual shall be forwarded to the subsequent owners or users of this device;
- (5) All specifications and design information in this document can be modified without notice. Some of the pictures in the user manual are schematic diagrams, which may not completely be consistent with the actual product.

1.2 Symbols

The descriptions of the symbols appearing in the interface, label, manual, and outer packaging of this device are as follows:

Symbols	Description	Symbols	Description
	Push-Pull Self-latching Connector (can be connected to an external foot switch)		Unsafe under Nuclear magnetic resonance examination environment
	SD Card interface		Fragile, handle with care
	USB3.0 interface		Keep dry
INPUT	Push-Pull Self-latching Connector		Temperature limit
	Power adapter interface		Date of manufacture
RESET	Reset button		Use-by date
VIDEO OUT	AV interface		Batch code
HDMI	HDMI interface		Serial number
	Refer to instruction manual/ booklet		Medical device
	Type BF applied part		The device should be sent to the special agencies according to local regulations
	Manufacturer		U.S. Federal law restricts this device to sale by or on the order of a licensed practitioner.

Symbols	Description	Symbols	Description
	Class II equipment		Caution
	General warning		Atmospheric pressure limitation
	Humidity limitation		

1.3 Equipment classification

- (1) Classification in accordance with the type of protection against electric shock: Class II equipment;
- (2) Classification in accordance with the degree of protection against electric shock: Type BF applied parts;
- (3) Classification in accordance with the degree of protection against liquid ingress: IPX0 ordinary equipment;
- (4) Classification in accordance with the degree of safety when used in the case of flammable anesthetic gas mixed with air or flammable anesthetic gas mixed with oxygen or nitrous oxide: Non-Type AP or APG equipment;
- (5) Classification in accordance with the operation mode: continuous operation;
- (6) Whether the equipment has an applied part to protect against defibrillation discharge effect: No;
- (7) Whether the equipment has a signal output or input portion: it has the signal output and input portions;
- (8) Permanently installed equipment or non-permanently installed equipment: non-permanently installed equipment;
- (9) Electromagnetic compatibility: Radiated emission and conducted emission are in line with Group 1 Class A of the CISPR 11 group classification requirements.

1.4 System safety requirements

- (1) Safety standards that the system meets:
 - IEC 60601-1 Medical electrical equipment Part 1: General requirements for basic safety and essential performance
 - IEC 60601-1-2 Medical electrical equipment Part 1-2: General requirements for basic safety and essential performance Collateral Standard: Electromagnetic disturbances - Requirements and tests
- (2) For detailed requirements of electromagnetic compatibility (EMC), please refer to "14 Attached Table".

2 Safety Guidelines

2.1 General warnings regarding use of device

- (1) This device must be used in accordance with this User Manual;
- (2) It must be used according to the intended purpose. Please refer to Section 3.4 Intended use;
- (3) Please prepare spare device that is immediately available for use and can be replaced in time in case of an emergency;
- (4) Please follow the User Manual of this device and associated device;
- (5) Keep the original packaging of the device in case you need to return it to the factory for repair;
- (6) It is forbidden to refit this device;
- (7) Always check the device before using it.
- (8) When using the device, it must be within your field of vision;

- (9) The operator cannot touch the patient and the device at the same time, and the patient cannot touch the device;
- (10) Improper use or maintenance of the device, exceeding the intended use, and premature wear or breakage of the device, can pose risks to patients and users;
- (11) This device must be used by trained personnel who have relevant knowledge and experience. Ensure that the user is able and ready to take appropriate action when the device encounters problems.

2.2 Warnings regarding installation and environment of use

- (1) This device is an MR unsafe device and cannot be used safely in an NMRI environment;
- (2) It is not recommended to disassemble the button battery, unauthorized disassembly may cause the risk of damage, if you need to disassemble the image processor, please contact Cook Medical;
- (3) When using this device, a proper space shall be reserved to remove the power plug from the side of the device to completely disconnect the device from the power supply;
- (4) Please turn off the device before disconnecting the supply mains.
- (5) The ME equipment could be isolated from supply mains by disconnection of the MAINS PLUG of the power supply cord of the AC Adaptor from supply mains.
- (6) Only operate this device under the specified environmental conditions;
- (7) Do not put anything on the cable. When laying out cables, be careful not to trip over others;
- (8) Do not bring this device near combustible materials;
- (9) The device shall be used to ensure adequate air circulation;
- (10) The device shall always be level during operation;
- (11) Do not use this device in the environment containing explosive or rich oxygen;
- (12) Make sure to use this device in a waterproof environment;
- (13) When transporting the device, packaging box shall be used to protect the device;
- (14) This device is not connected to the ultrasonic endoscope and is not the power supply device thereof;
- (15) This device is not connected with the endoscopic accessories of extracorporeal lithotripsy and their medical electrical equipment;
- (16) This device is not connected with the tissue's ultrasound-attracted endoscopic accessories and their medical electrical equipment.

2.3 Warnings regarding electromagnetic interference

- (1) Ensure that all equipment operating in the vicinity meet the relevant electromagnetic compatibility (EMC) requirements. Electromagnetic radiation from connected external equipment (such as video equipment) can affect image quality. In the case of very serious electromagnetic interferences, the image quality will be severely disturbed (for example, slight streaks and color changes appear on the display);
- (2) The Endoscopic Video Image Processor complies with the electromagnetic compatibility requirements in the IEC 60601-1-2 standard, and the user shall install and use it according to the electromagnetic compatibility information provided;
- (3) Medical electrical equipment is subject to specified precautions for electromagnetic compatibility (EMC). Please follow the Precautions on Electromagnetic Compatibility (EMC) of Medical Electrical Equipment

provided with the User Manual;

- (4) Portable and mobile RF communication equipment may affect the performance of Endoscopic Video Image Processor. Avoid strong electromagnetic interference when using, such as near mobile phones, microwave ovens, etc;
- (5) Class A equipment is intended for use in industrial environments. Due to conduction disturbance and radiation disturbance from the Endoscopic Video Image Processor, it may be potentially difficult to ensure electromagnetic compatibility in other environments;
- (6) Equipment or systems shall not be used close to or stacked with other equipment. If it must be used close to or stacked with other equipment, it shall observe and verify that it can work under the configuration used;
- (7) Leakage current multiplication - patient risk: In case that the product is used in combination with medical electrical devices and/or motorized endoscopic accessories, the leakage current may be multiplied;
- (8) If this device or endoscopic accessories are used with devices from different manufacturers, as well as with medical electrical devices, it shall be ensured that they meet the BF conditions (insulated, ungrounded applied part), that is, devices with the identification; 
- (9) Don't near active HF surgical equipment and the RF shielded room of an ME system for magnetic resonance imaging, where the intensity of EM disturbances is high;
- (10) Use of accessories, transducers and cables other than those specified or provided by the manufacturer of this equipment could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation;
- (11) Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the equipment, including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result;
- (12) The emissions characteristics of this equipment make it suitable for use in industrial areas and hospitals (CISPR 11 class A). If it is used in a residential environment (for which CISPR 11 class B is normally required), this equipment might not offer adequate protection to radio-frequency communication services. The user might need to take mitigation measures, such as relocating or re-orienting the equipment.
- (13) Electromagnetic interference* (EMI) may occur on this instrument near equipment marked with the following

symbol  . Electromagnetic interference may occur on this instrument near other portable and mobile radio frequency (RF) communications equipment or specific RF emitters for example, cellular phones, diathermy, electrocautery, and RFID) or other specific potential electromagnetic disturbance sources. To check for EMI, verify the system's operation in which it will be used. Should EMI occur, employ mitigation measures like reorienting or repositioning the instrument, or shielding its location. Placing this instrument near other medical electrical equipment or mobile RF communications equipment may result in EMI, which may degrade the video image.

2.4 Warnings regarding cable and accessories

- (1) The device and/or accessories shall be checked before use and ensure that they are functioning properly;
- (2) It must be verified that the accessories and/or peripheral equipment connected to the interface of this device comply with the relevant specifications (IEC 60601-1). Moreover, all the configurations of the system must

- meet the requirements of IEC 60601-1-2;
- (3) Failure of operation may result from damage to the internal electrical components of the device due to pulling of video cables or heavy pressure on the device;
 - (4) All electrical connections shall be checked regularly. Loose plugs, video cable failure, Endoscopic Video Image Processor damage and other issues will affect the image quality and shall be replaced immediately;
 - (5) Hold the plug and pull it out of the power supply. Do not pull the cable;
 - (6) Only use the cable provided with the equipment to operate this device;
 - (7) When handling the patient do not simultaneously touch the power socket or SIP/SOP connector of Endoscopic Video Image Processor;
 - (8) The Push-Pull Self-latching Connector of the Endoscopic Video Image Processor is only applicable to the endoscope distributed by Cook Medical, and cannot be mixed with other manufacturers;
 - (9) Do not fire the laser until the laser probe is put in the correct position displayed in the endoscopic image, keep a proper distance between the target and the tip of the endoscope, and use the lowest possible power output in case of patient injury, burn, bleeding, perforation, or device damage.

3 Device Description

3.1 Name

Endoscopic Video Image Processor

3.2 Model

ASUS-IPD-V100

3.3 Structure

This device is composed of image processor (including software), power cord, power adapter, video cable (image transmission) , HDMI cable(optional) and foot switch (optional).

3.4 Intended use

The device has been developed to be used with Single-Use Flexible Ureteroscope and endoscopic accessories such as a biopsy forceps for endoscopy and endoscopic surgery within urinary tract and interior of the ureter and kidney.

3.5 Contraindications

This device is a supporting device of the video endoscopes distributed by Cook Medical. For related contraindications, please refer to the Ascend Single-Use Flexible Ureteroscope.

4 Control Operation of the Device Panel

Note: All pictures shown are for reference only, please refer to the actual object. For all the pictures mentioned below, please refer to the page 1 of the manual.

4.1 Front view of the Endoscopic Video Image Processor (Figure 1)

4.2 Rear view of the Endoscopic Video Image Processor (Figure 2)

4.3 Installation guide(as shown in the Figure 3)

Step 1: Open the DC&RESET interface cover on the left side of the Endoscopic Video Image Processor.

Step 2: Insert the DC connector of the power adapter into the DC interface.

Step 3: Insert the “8”-shaped connector of the power cord into the socket of the power adapter.

Step 4: Install the video endoscope.

- (1) Align the interface of the video cable with the interface of handle;
- (2) Insert the video cable into the handle interface;
- (3) Align the insertion mark of the Endoscopic Video Image Processor connector;
- (4) Insert the Endoscopic Video Image Processor connector;

(5) The video cable connection is completed.

Step 5: After installation, check that there is no plug missing, and then plug the other end of the power cord into the power socket to make sure it is connected to the power supply.

The original video cable, power cord and power adapter provided by Cook Medical must be used.

5 Basic Operation

5.1 Startup and Login

Press the Power button  of the Endoscopic Video Image Processor and select the appropriate language (as shown in the Figure 4). Boot login interface (as shown in the Figure 5) factory default condition. User name: admin, password: 888888.

After entering the default password, the system will enter into administrator mode. If you choose “Skip”, you will skip this login interface and enter main interface in guest mode. If you select “This page will not be displayed after startup”, the default is entering main interface when it is turned on next time and no password is required to log in again.

The function permissions of the two modes are shown in the table below:

Function Access	Administrator Mode	Guest Mode
Patient information	All functions	Query, create new patient information, query patient images
Endoscopy	All functions	All functions
Menu Settings	All functions	Functions other than user management, image settings (except for frame shape adjustment), factory reset, software version number, and software upgrade.

Note:

- (1) If you check the option "This page will not be displayed after startup", you will no longer need to enter a password for subsequent startups.
- (2) If you want to perform a function that is not authorized in the guest mode, you will be prompted to enter the administrator password. After the password is entered, it will automatically switch to the administrator mode. In the main interface of real-time video, click the floating ball to enter the main interface. (as shown in the Figure 1)

5.2 Real-time video

- (1) Click the **Video**  button from the main interface to enter the Real-time video interface. (as shown in the Figure 6)
- (2) Click the button from the main interface of real-time video to create new patient information. Click the button to save the information. (as shown in the Figure 7)
- (3) Click the button from the main interface of real-time video to query patient information. Click the button to search for patient information. (as shown in the Figure 8)
- (4) Click the corresponding information column of the patient you want to view and enter the detailed information interface of the patient. (as shown in the Figure 9)

- (5) Click the  button in the upper right corner to apply the corresponding patient information and return to the Real-time video interface.(as shown in the Figure 9)
- (6) Click the  button in the upper right corner to delete the patient information.(as shown in the Figure 9)
- (7) Click the  button in the upper right corner to export the patient file. For image/video playback, see 5.7. (as shown in the Figure 9)

5.3 White balance

In the main setting interface(as shown in the Figure 12), you can set “Fix WB”(Fixed White Balance) and “Manual WB”(Manual White Balance). The fixed white balance has been adjusted in the factory settings, and users are not allowed to debug. The manual white balance can be adjusted by the user by pressing the **WHT BAL**  button.

The manual white balancing steps are as follows:

- ① Aim the endoscope vertically at a white object, such as white gauze, white paper.
- ② Press the **WHT BAL**  button and keep aiming at the white object, the white balance can be completed after 2~3 seconds.

5.4 Adjust the screen brightness/image display area size

- ① Adjust screen brightness

To adjust the screen brightness, press the **BRIGHTNESS +**  / **BRIGHTNESS -**  button until the desired brightness is reached. There are 8 brightness levels to adjust.

- ② Adjust the size of the image display area

Adjust the size of the image display area (default: 800×800), press the **ZOOM in/out**  button (4 levels in total), and select the best image observation size.

5.5 Take a picture/freeze the picture

- ① In the real-time video interface with applying patient information, press the **PHOTO**  button to freeze the screen. The prompt “Screen Freeze” appears on the top of the screen which means it is in the frozen state.
- ② When it is in the frozen state, press the **PHOTO**  button again and the screen freeze will be released.
- ③ When it is in the frozen state, long press the **PHOTO**  button for 3 seconds to save the frozen picture. Then the prompt "Save Successfully" appears at the top of the screen, the picture is taken. If device connects to the foot switch, short press the left switch also can take a picture.
- ④ The maximum resolution of the stored images: 800*800. The images can be played back from patient information. **For image viewing/video playback, see 5.7.**

Do not use the image processor if real-time endoscopic video cannot be observed. Otherwise, the patient may be injured.

5.6 Record

- ① In the real-time video interface with applying patient information, press the **VIDEO**  button to start recording. The prompt “Recording XX:XX:XX” appears on the top of the screen which means it is in the record state. (as shown in the Figure 10)
- ② When it is in the record state, press the **VIDEO**  button again to complete the recording and save the video.
- ③ If device connects to the foot switch, short press the right switch to start recording and press again to save the video and exit recording.
- ④ The maximum resolution of the stored videos: 800*800. The videos can be played back from patient information. **For image viewing/video playback, see 5.7.**

Note: During recording, the image cannot be rotated by clicking the screen or pressing the foot switch.

5.7 Image viewing/video playback

- ① Click Query to enter the patient information interface (this interface can be seen in 5.2), select a row of target patient information, enter the detailed patient information interface (this interface can be seen in 5.2), select a picture, and the image can be played back (as shown in the Figure 11)
- ② Select the video and click the **Play** icon to play the video. (as shown in the Figure 11)

5.8 Menu Setting

Click the **Menu**  button from main interface to enter the menu setting main interface. (as shown in the Figure 12)

- ① Language Settings: Click the language button to select your preferred language.
- ② User Management: Administrator can enter the user management interface to perform password operations by clicking the button.
- ③ Image Settings: There are 5 modes of image setting: Custom 1, Choleidoscope, Hysteroscope, Cystoscope and Ureteroscope. The levels of image parameters setting are shown as numbers (as shown in the Figure 13).

After saving preferred image settings, return to the Real-time video interface, press the **USER MODE**  button of device panel can apply corresponding mode.

- ④ Time Settings: Click the button to set the current date and time.
- ⑤ Resolution: The adjustment is the initial size of the image when entering the Real-time video interface. There are 2 resolutions available for selection: 400*400 and 800*800.
- ⑥ Factory Reset: Administrator can click the button to reset the system to factory defaults.
- ⑦ Software Version: Review the current software version. Continuously click the Software Version for 5 times to enter the Engineering Mode.
- ⑧ Software Update: Please contact Cook Medical if a software update is required.


Warning: This function is limited to authorized personnel.

- ⑨ Image Records of patient Information: Click the  button, and the patient information will be displayed on the picture when taking photos.

5.9 Image rotation

When observing an image or taking a picture, you can click the **Angle Rotation** button at the bottom right corner of the screen to rotate the image clockwise (cycle among 0°, 90°, 180° and 270°), as shown in the figure below.

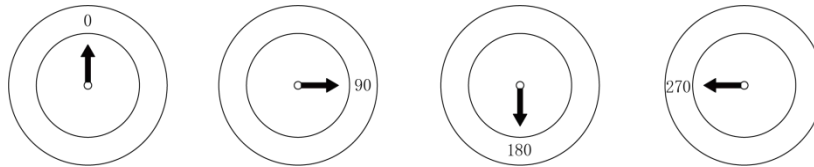Note: To prevent injury to the patient, use this function with caution.

5.10 Turn off

When turning off the device, the user can choose whether to exit the administrator mode. (as shown in the Figure 14)

5.11 Foot switch (optional) function

There is a left switch and a right switch on the foot switch. (as shown in the Figure 15)

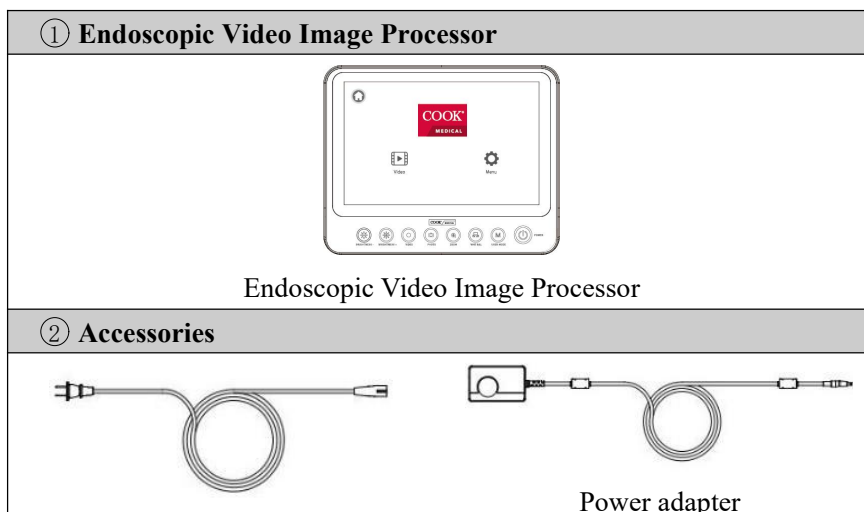
Right switch: short press to start recording, press again to save the video and exit recording.

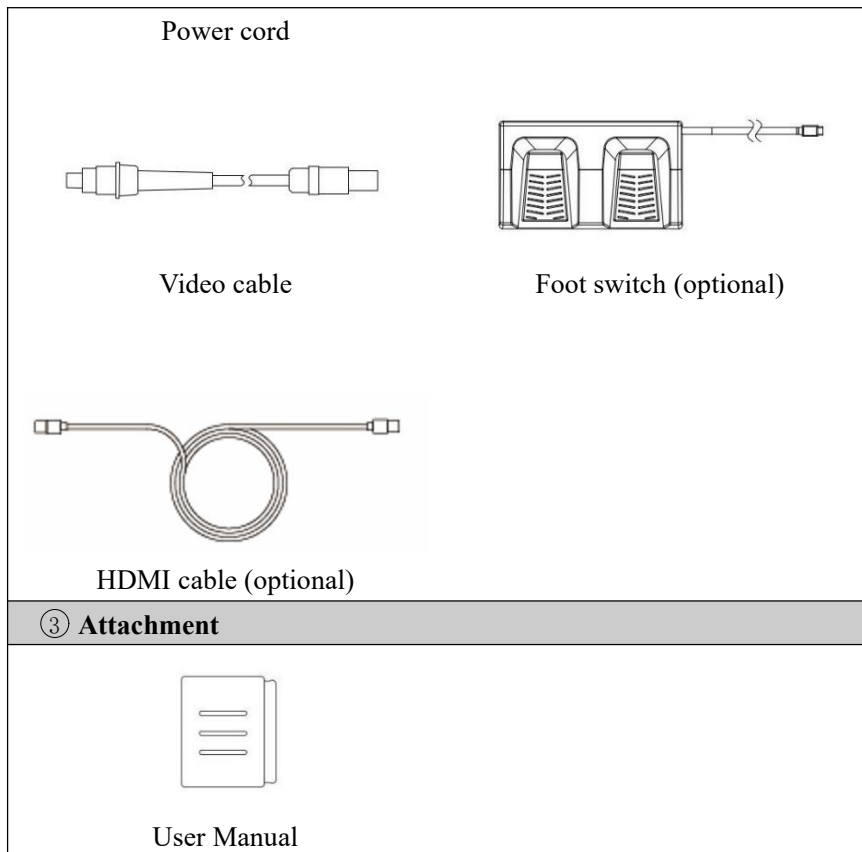
Left switch: short press to take a picture. Long press to rotate the image clockwise. (Cyclic rotation at 0°, 90°, 180° and 270°)

That before each use or after a change of viewing modes/settings, the OPERATOR should check to ensure the view observed on screen through the ENDOSCOPE provides a live image (rather than a stored one) and has the correct image orientation.

6 List of Supplies

7 Precautions for Use

If the user fails to comply with the following guidelines, the manufacturer will not be held responsible for the possible consequences.

Before putting the equipment into use, please make sure that the electrical connection complies with:

- (1) Appropriate technical regulations; appropriate fire and explosion safety regulations.
- (2) The safety of users and patients depends on whether the cable is in good condition.
- (3) The cable must be checked before use, otherwise it must not be used.

Other equipment required for assembly/initial installation:

Our Single-use flexible ureteroscopes is a single-use sterile device. Before using it for the first time, make sure that the packaging is intact. If the packaging is damaged, it is forbidden to use it. After use, please dispose of this device in accordance with the hospital's waste disposal procedures, and it is forbidden to reuse it.

Make sure that the video cable is not broken, knotted or twisted. If the video cable is damaged, do not continue to use it.

Make sure that the pins in the plug are not bent and that the plug is not wet.

Cable

Check whether the cables are in good condition and that there are no cracks, knots or frays.

USB

The USB flash disk to be used for interaction must be free of viruses and Trojan horses.

External Monitor

Before connecting the external monitor, please ensure that the external monitor have HDMI or DVI input signal type; and it supports 1080P resolution or it has adaptive function: the external monitor can receive 1080P images internally and convert them to physical resolution. In addition, the recommended specifications of external monitor are as follows:

- luminance: $\geq 250\text{cd/m}^2$
- color gamut: $\geq 72\%\text{NTSC}$
- contrast ratio: $\geq 500:1$
- monitor size: ≥ 17 inch

8 Cleaning, Disinfection and Maintenance

You should use cleaning agents and disinfectants that are approved for use in this device. Follow the recommended concentration, temperature, action time and use time as below:

- (1) Only wipe the device with a damp cloth. Do not immerse the device in liquid.
- (2) Do not use abrasive cleaners or solvents, otherwise the surface or label will be damaged.
- (3) Do not use wet sponges or wet cloths. Too much cleaning agent will come into contact with electronic components and cause equipment damage.
- (4) Only after the cleaned components are completely dry, can the system be reconnected to the power supply.
- (5) Endoscopic Video Image Processor is a non-sterilized device. Do not soak or sterilize the processor, otherwise it will be damaged.

Note:

The device needs to be cleaned, inspected and maintained every other week as below:

- (1) Cleaning: Please wipe the outer surface of the instrument (including Power cord, Power adapter, Endoscopic Video Image Processor shell, Endoscopic Video Image Processor screen) with a dust-free cloth.
- (2) Disinfection: After the instrument is cleaned, wipe the outer surface of the instrument with 3-4 drops of 75% alcohol with a dust-free cloth.
- (3) If you encounter stubborn stains, wipe it with 1-2 drops of 75% alcohol with a dust-free cloth, until there is no stain on the surface.
- (4) The device shall be cleaned after each use.
 - ① Power cord: Before and after each use, check whether the cord is damaged or broken.
 - ② Power adapter: Before and after each use, check for damage or breakage.
 - ③ Endoscopic Video Image Processor shell: Before and after each use, check whether the shell is damaged.
 - ④ Endoscopic Video Image Processor screen: Before and after each use, please wipe the surface of the instrument with a dust-free cloth to check for damage.

9 Troubleshooting

If the processor is visibly damaged, does not function as expected, do not use the processor.

Some problems that appear to be malfunctions may be correctable by referring to the section. If the problem cannot be resolved by the described remedial action, stop using the processor.

Fault	Possible cause	Solution
The image is too dark and the illumination brightness is too low	Not adjusted to the appropriate level of brightness	Press the Brightness +/- button to adjust to the appropriate brightness level.
No image on the display	Not powered	Turn on the Endoscopic Video Image Processor.
	The power adapter interface is loose or incorrectly connected	Connect the instrument according to Chapter 4.3 Installation Guide, and then turn on the power switch.
Color bars are displayed on the display	The video cable is not plugged in	Insert the video cable correctly.
	Video cable failure	Replace the video cable or send it for repair.
Unnatural color	The white balance is not performed correctly.	Perform white balance again, see the Chapter Basic Operation.
	The color of the display is not adjusted properly reset.	Restore factory settings.
	The video cable is faulty	Replace the video cable.
Image distortion	Incorrect settings	Follow the Chapter Assembly/Initial Installation for correct settings.

10 Guarantee

- (1) The manufacturer guarantees the quality of the device for 12 months. This guarantee is limited to claims made in written form from the invoice date and within the warranty period, subject to the invoice number. The legal guarantee statement is no longer within the limits of this guarantee;
- (2) This guarantee does not apply to failures caused by daily wear and tear, improper use, improper operation, poor maintenance or natural forces;
- (3) If the device is repaired by the user, the guarantee and warranty statement will be invalid. If the equipment needs maintenance, the above requirements also apply;
- (4) No claims can be made for problems caused by incorrect operation or use with other equipment or accessories.

11 Discard

- (1) Some risks of infection might be caused by the contaminated device, hence it must be washed and disinfected before disposal;
- (2) When discarding or reusing this device and its accessories, it should be disposed of in accordance with local laws or medical waste;
- (3) The cleaned and disinfected device should be stored in a package that can isolate bacteria, where is placed in a dark, cool, and appropriate room;
- (4) When used in accordance with the provisions of the user manual, the service life of this product is 5 years.

12 Basic Parameters

12.1 Basic parameters

Item name		Description
Product model		ASUS-IPD-V100
Operating system		Linux 4.9.65
Memory		Micro SD 32G
Power adapter		12V 5A
Voltage		AC 100V-240V
Screen parameters	Screen size	10.1" LCD
	Resolution	1280 × 800 pixels
	Touch screen	Capacitive touch screen, and multi-point touch
Interface parameters¹⁾	DC	12V 5A
	HDMI	Connect to the display to observe the image.
	USB	Connect the USB to export the file, the format of the exported file is video (MP4 format) and picture (BMP format).
	Push-Pull Self-latching Connector 2	Connect the endoscope to achieve foot switch control.
	AV (CVBS)	Connect to the display to observe the image
	SD card	Insert the SD Card into the Endoscopic Video Image Processor to store the image. The information in the SD Card can be read on the computer.
	Push-Pull Self-latching Connector	Connect the video endoscope for image signal transmission
Device size		264mm*210mm*38mm
Video cable (image transmission)		Connector at one end, Push-Pull Self-latching Connector at the other end
Software version		V2
Adaptable devices		Suitable for video endoscopes distributed by Cook Medical

¹⁾ Additional equipment connected to medical electrical equipment through the interfaces must comply with the respective IEC or ISO standards (e.g. IEC 60950 or IEC 62368 for data processing equipment). Furthermore all configurations shall comply with the requirements for medical electrical systems (see clause 16 of the latest valid version of IEC 60601-1). Anybody connecting additional equipment to medical electrical equipment configures a medical system and is therefore responsible that the system complies with the requirements for medical electrical systems. If in doubt, consult your local representative or Cook Medical.

12.2 Environmental requirements

Item	Operating environment	Transport environment	Storage environment
Ambient temperature	10°C - 40°C	-20°C - 55°C	10°C - 40°C
Relative humidity	30%-85%	30%-85%	30%-85%
Atmospheric pressure	70kPa-106kPa	70kPa-106kPa	70kPa -106kPa

13 Software-related Instructions

13.1 Operating environment

13.1.1 Hardware configuration

Processor: TI AM5708

Storage: SD

13.1.2 Software environment

System software: Linux

13.1.3 Network conditions

Network card: None

Network type: None

Network architecture: None

13.2 Security software: No security software.

13.3 Data and equipment (system) interface

Transmission protocol: HDMI interface (TMDS minimized transmission differential signal protocol), AV interface (composite video transmission protocol), USB interface (Universal Serial Bus), Push-Pull Self-latching Connector (Analog Signal), and SD Card interface (SDIO Secure Digital Input and Output Protocol)

13.4 User access control mechanism

User authentication method: After booting, the interface is the user login interface. Login can choose login and administrator login according to user name and password.

User type and authority: The system includes two user levels: administrator and ordinary user. The administrator can modify user information after logging in; ordinary users can operate the functions displayed by the operating system.

13.5 Software environment (including system software, supporting software, application software) and relevant requirements for security software updates.

The software updates are only to be performed by trained and authorized personnel.

14 Attached Table

Table 1 Guidance and manufacturer's declaration - electromagnetic emissions

Emissions test	Compliance
RF emissions CISPR 11	Group 1
RF emissions CISPR 11	Class A
Harmonic emissions IEC 61000-3-2	Not Applicable
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Not Applicable

Table 2 Guidance and manufacturer's declaration - electromagnetic immunity

Immunity Test	IEC 60601-1-2 Test level	Compliance level
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air
Electrical fast transient/burst IEC 61000-4-4	±2 kV power supply lines ±1 kV signal input/output 100 kHz repetition frequency	±2 kV power supply lines ±1 kV signal input/output 100 kHz repetition frequency
Surge IEC 61000-4-5	±0.5 kV, ±1 kV differential mode±0.5 kV, ±1 kV, ±2 kV common mode	±0.5 kV, ±1 kV differential mode ±0.5 kV, ±1 kV, ±2 kV common
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0 % U _T ; 0,5 cycle. At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315°. 0 % U _T ; 1 cycle and 70 % U _T ; 25/30 cycles; Single phase: at 0°. 0 % U _T ; 250/300 cycle	0 % U _T ; 0,5 cycle. At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315°. 0 % U _T ; 1 cycle and 70 % U _T ; 25/30 cycles; Single phase: at 0°. 0 % U _T ; 250/300 cycle
Power frequency magnetic field IEC 61000-4-8	30 A/m 50Hz/60Hz	30 A/m 50Hz/60Hz
Conducted RF IEC61000-4-6	3 V 0,15 MHz – 80 MHz 6 V in ISM bands between 0,15 MHz and 80 MHz 80 % AM at 1 kHz	3 V 0,15 MHz – 80 MHz 6 V in ISM bands between 0,15 MHz and 80 MHz 80 % AM at 1 kHz
Radiated RF IEC61000-4-3	3 V/m 80 MHz – 2,7 GHz 80 % AM at 1 kHz	3 V/m 80 MHz – 2,7 GHz 80 % AM at 1 kHz
NOTE U _T is the a.c. mains voltage prior to application of the test level.		

Table 3 Guidance and manufacturer's declaration - electromagnetic immunity

	Test Frequency (MHz)	Band (MHz)	Service	Modulation	IEC 60601-1-2 Test Level (V/m)	Compliance level (V/m)
Radiated RF IEC61000-4-3 (Test specifications for ENCLOSURE PORT IMMUNITY to RF wireless communications equipment)	385	380 –390	TETRA 400	Pulse modulation 18 Hz	27	27
	450	430 –470	GMRS 460, FRS 460	FM ± 5 kHz deviation 1 kHz sine	28	28
	710	704 – 787	LTE Band 13, 17	Pulse modulation	9	9
	745					

	780	800 – 960	GSM 800/900, TETRA 800, iDEN 820, CDMA 850, LTE Band 5	217 Hz	28	28
	810			Pulse modulation 18 Hz		
	870					
	930					
	1720	1700 – 1990	GSM 1800; CDMA 1900; GSM 1900; DECT; LTE Band 1, 3, 4, 25; UMTS	Pulse modulation 217 Hz	28	28
	1845					
	1970					
	2450	2400 – 2570	Bluetooth, WLAN, 802.11 b/g/n, RFID 2450, LTE Band 7	Pulse modulation 217 Hz	28	28
	5240	5100 – 5800	WLAN 802.11 a/n	Pulse modulation 217 Hz	9	9
	5500					
	5785					

Table 4 Guidance and manufacturer’s declaration - electromagnetic immunity

Radiated RF IEC61000-4-39 (Test specifications for ENCLOSURE PORT IMMUNITY to proximity magnetic fields)	Test Frequency	Modulation	IEC 60601-1-2 Test Level(A/m)	Compliance level (A/m)
	30 kHz	CW	8	8
	134,2 kHz	Pulse modulation 2.1 kHz	65	65
	13,56 MHz	Pulse modulation 50 kHz	7,5	7,5

- (1) Except for the adapters, power cords, and video cables sold by the manufacturer of the equipment or system as spare parts for internal components, the use of accessories, adapters, and video cables outside the regulations may lead to an increase in the emission of the equipment or system or a decrease in the immunity. The following cable types must be used to ensure compliance with interference radiation and immunity standards:

Cable	Video cable	Power adapter cable	Power cord	Foot switch wire
Length (m)	2.8	1.5	2.0	1.8

- (2) Please use the following power adapters:
Manufacturer: BRIDGEPOER CORP, Model: MENB1060A1203N01

- (3) You should use the HDMI cable:
Manufacturer: UGREEN GROUP LTD, Model: HD156
- (4) The HDMI interface of the device can be connected with the HDMI cable externally and connected to the display (not part of the device structure) for image observation. In order to meet the standards of interference radiation and immunity, the HDMI cable with the length of less than 3m is required
- (5) When connecting the external display with HDMI interface, you should use the HDMI cable recommended by the manufacturer.
- (6) Use of this equipment adjacent to or stacked with other equipment should be avoided because it could result in improper operation. If such use is necessary, this equipment and the other equipment should be observed to verify that they are operating normally.
- (7) Active medical devices are subject to special EMC precautions, and therefore must be installed and used in accordance with these guidelines.
- (8) Portable and mobile communication radio frequency equipment may affect the use of medical electrical equipment.
- (9) Basic performance description: Before, during and after the test, check the imaging function. The device can work and the noise in the image shall not affect the diagnosis conclusion; during the test, the Endoscopic Video Image Processor shall be able to maintain its set working status, without crashes, working mode changes or other malfunctions, the white balance, the camera and video functions are normal.
- (10) Even if other equipments meet the emission requirements of corresponding national standards, the equipment or system may be interfered by other equipment.

Manufactured By:

Guangzhou Red Pine Medical Instrument Co.,Ltd.

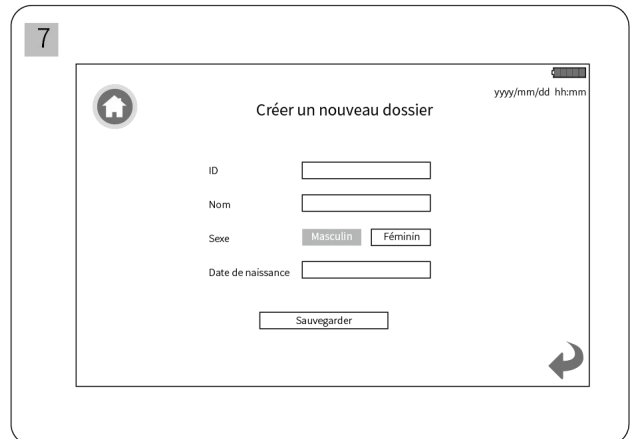
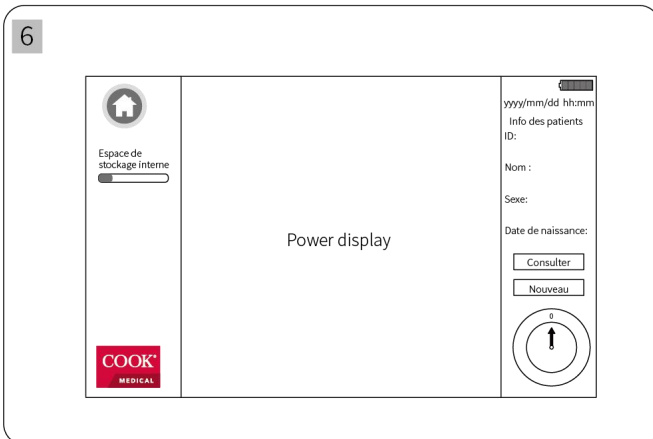
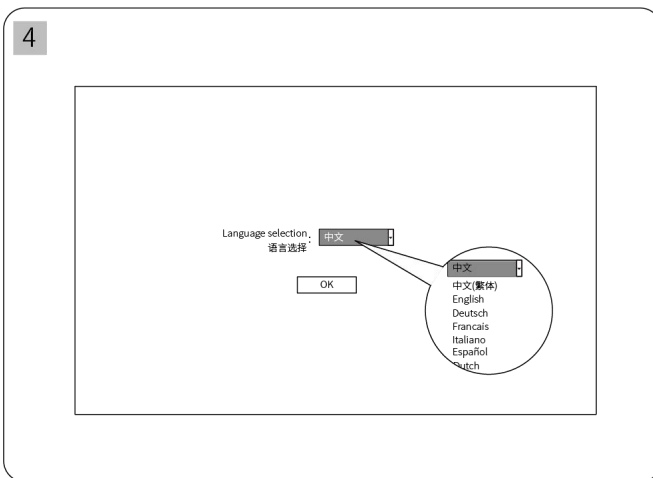
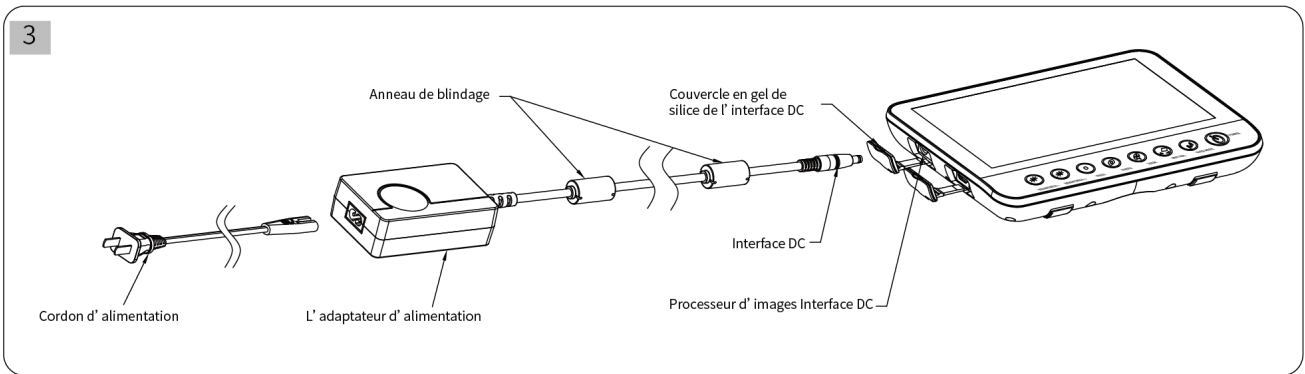
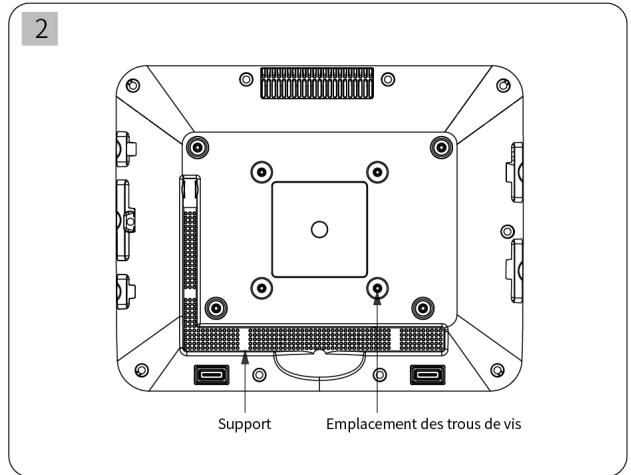
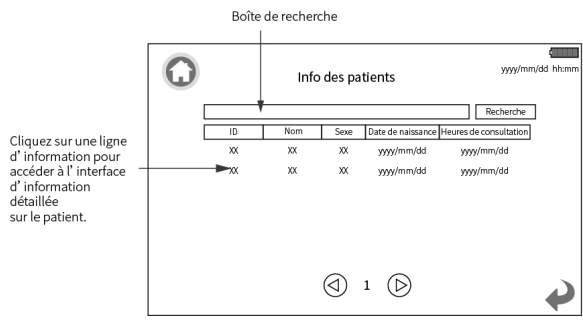
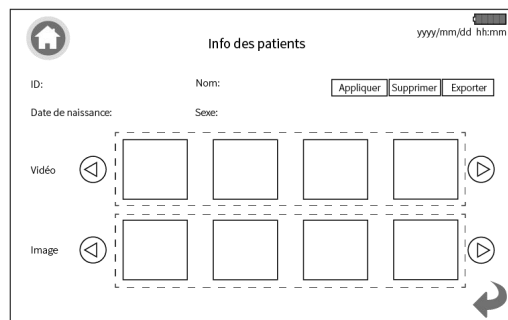
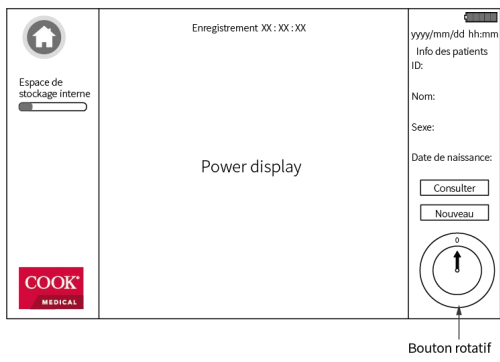
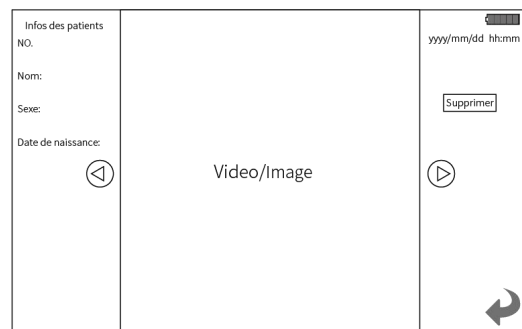
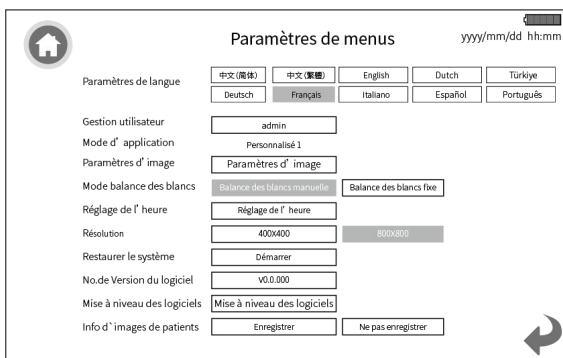
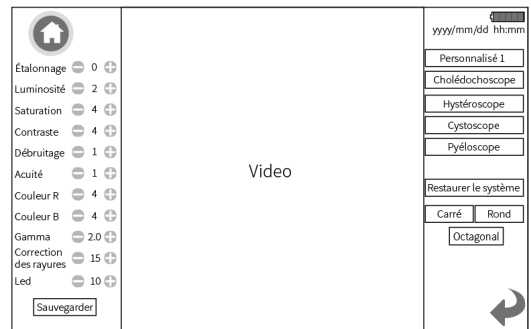
12 F, No.87 Luoxuan Avenue, Guangzhou International Bioisland, Huangpu District, 510000 Guangzhou City,
Guangdong Province, PEOPLE'S REPUBLIC OF CHINA

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Fabricant: UGREEN GROUP LTD, Modèle: HD156
Guangzhou Red Pine Medical Instrument Co., Ltd.
WILLIAM A. COOK AUSTRALIA PTY. LTD.