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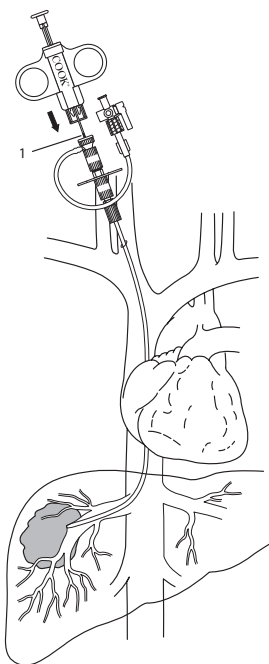
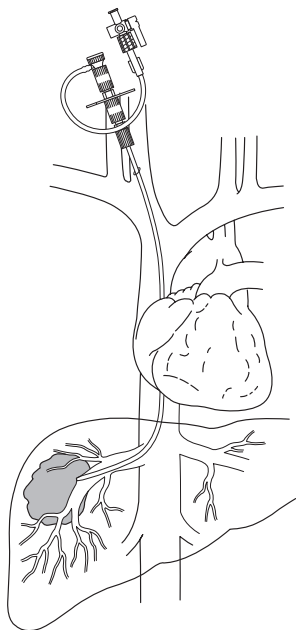
Fig. 3

LIVER ACCESS AND BIOPSY SET

CAUTION: U.S. federal law restricts this device to sale by or on the order of a physician (or properly licensed practitioner).

DEVICE DESCRIPTION

The Liver Access and Biopsy Set includes a Quick-Core® Biopsy Needle (**Fig.1**) and a preassembled liver access set consisting of a radiopaque sheath, a Check-Flo® valve adapter and a stiffening cannula. The set also may contain one straight catheter, one curved catheter, and/or one dilator.

INTENDED USE

The Liver Access and Biopsy Set is intended for use in obtaining liver histology samples via a jugular vein approach. The transjugular approach provides an alternative to large-needle percutaneous biopsy in patients with abnormal clotting factors or in patients with ascites, which may decrease the success of routine percutaneous biopsy and increase the risk of hemorrhagic complications.

The Quick-Core Biopsy Needle is intended for soft tissue biopsy.

CONTRAINDICATIONS

None known

WARNINGS

Extreme care must be exercised during manipulation and withdrawal of catheter to prevent pulling catheter apart. Insertion through a synthetic vascular graft should be avoided whenever possible.

PRECAUTIONS

- The product is intended for use by physicians trained and experienced in transjugular liver biopsy techniques. Standard techniques for transjugular liver biopsy should be employed.
- Fluoroscopic or ultrasound visualization is required to ensure accurate placement.
- **When test firing Quick-Core Biopsy Needle, place distal beveled tip of inner stylet against solid surface before activating cutting cannula. This will prevent damage to cutting cannula tip.**
- The potential effects of phthalates on pregnant/nursing women or children have not been fully characterized and there may be concern for reproductive and developmental effects.

POTENTIAL ADVERSE EVENTS

- Intraperitoneal hemorrhage
- Puncture site hematoma
- Cardiac arrhythmias
- Arteriovenous fistulas

INSTRUCTIONS FOR USE

1. Utilizing standard access techniques, introduce an appropriate .035 inch wire guide into the inferior vena cava via jugular vein access.
2. Using a selective catheter and wire guide of choice, catheterize the right hepatic vein or an appropriate alternative hepatic vein branch. Leave the wire in a safe, distal position, and remove the catheter. **CAUTION: To prevent cardiac arrhythmias, continuous cardiac monitoring is recommended while negotiating past the right atrium.**
3. Introduce and advance the transjugular introducer sheath and stiffening cannula over the wire guide into the selected hepatic vein. When introducing these components as a preassembled set, the included straight catheter may be used to facilitate introduction. Once the set is positioned within the hepatic vein, the straight catheter should be removed. **NOTE:** Care should be taken to prevent damage to the straight catheter during removal through the metal stiffening cannula. Leaving a wire guide through the straight catheter during removal may aid in preventing catheter damage.
4. Prepare the Quick-Core Biopsy Needle by pulling back on the plunger until a firm click is felt, indicating the needle spring is locked into the ready position.
5. With the stylet fully retracted, so that the specimen notch is completely covered by the cutting cannula, advance the Quick-Core Biopsy Needle through the liver access set until the etch mark on the proximal cannula is at the hub of the Check-Flo valve adapter. (**Fig. 2**) This will position the needle tip at the distal end of the access assembly.
6. Using the directional arrow on the hub of the stiffening cannula, direct the liver access set in the desired direction for biopsy and apply gentle forward pressure to "tent" the wall of the hepatic vein adjacent to the biopsy site. (**Fig. 3**)
7. Advance the Quick-Core Biopsy Needle out the tip of the liver access assembly and into the liver tissue. **Do not advance the stylet until the needle is in position.**
8. While maintaining needle position, advance the stylet to expose the specimen notch within the area to be biopsied. Fire the cutting cannula

by fully depressing the plunger to capture tissue within the specimen notch.

9. Withdraw the Quick-Core Biopsy Needle from the biopsy area and completely out of the liver access set.
10. To remove the tissue specimen, pull back on the plunger until a firm click is felt. This indicates that the cutting cannula is locked into position. Push the stylet forward to expose the tissue specimen and remove the specimen from the notch.

HOW SUPPLIED

Supplied sterilized by ethylene oxide gas in peel-open packages. Intended for one-time use. Sterile if package is unopened or undamaged. Do not use the product if there is doubt as to whether the product is sterile. Store in a dark, dry, cool place. Avoid extended exposure to light. Upon removal from package, inspect the product to ensure no damage has occurred.

REFERENCES

These instructions for use are based on experience from physicians and (or) their published literature. Refer to your local Cook sales representative for information on available literature.

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© COOK 2019
Bloomington, IN 47404 U.S.A.
Cook Medical Europe Ltd