COOK MEDICAL

CE 0123

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ENGLISH

TABLE OF CONTENTS

1. DEVICE DESCRIPTION 23
- 1.1 Construction of Device 23
- 1.2 Performance Characteristics 23
- 1.3 Specifications 23
- 1.4 Device Compatibility 24
- 1.5 Qualitative and Quantitative Information 24
- 1.6 Patient Population 24
- 1.7 Intended User 24
- 1.8 Contact with Body Tissue 24
- 1.9 Operating Principle 24
- 1.10 Expected Lifetime of the Device 24

2. INTENDED USE 25

3. INDICATIONS FOR USE 25

4. CLINICAL BENEFITS 25

5. CONTRAINDICATIONS 25
- 5.1 Filter Placement 25
- 5.2 Optional Filter Retrieval 25

6. WARNINGS 25
- 6.1 Sterile and Single Use 25
- 6.2 Malfunctions/Changes in Performance 25
- 6.3 Optional Filter Retrieval 26

7. PRECAUTIONS 26
- 7.1 General Precautions 26
- 7.2 Malfunctions/Changes in Performance 26
- 7.2.1 Malfunctions/Changes in Performance for Femoral Filter Placement 26
- 7.2.2 Malfunctions/Changes in Performance for Jugular Filter Placement 26
- 7.3 CMR 1A/1B and/or Endocrine Disrupting Substances 27
- 7.4 Allergic Reaction 27
- 7.5 Training Requirements for User 27
- 7.6 Requirements for Facilities 27

8. MRI SAFETY INFORMATION 28

9. POTENTIAL ADVERSE EVENTS 29

10. STERILIZATION 29

11. HOW SUPPLIED 29

12. INSPECTION OF DEVICE 29

13. DEVICE PREPARATION 29
- 13.1 Device Selection 29

14. INSTRUCTIONS FOR USE 29
- 14.1 Patient Preparation 29
- 14.2 Step-by-Step Guide 29
- 14.2.1 Femoral Approach 29
- 14.2.1.1 Preparation 29
- 14.2.1.2 Filter Placement 29
- 14.2.2 Jugular Approach 30
- 14.2.2.1 Preparation 30
- 14.2.2.2 Filter Placement 30
- 14.2.3 Optional Retrieval Procedure 31

15. REQUIRED EQUIPMENT AND ACCESSORIES 31
- 15.1 Materials Required 31
- 15.2 Materials Recommended 31

16. DISPOSAL OF DEVICE 31

17. REFERENCES 32

18. PATIENT COUNSELING 32
- 18.1 Symptoms 32
- 18.2 Follow Up Information for the Patient 32
- 18.3 Patient Implant Information 32
- 18.4 Summary of Safety and Clinical Performance (SSCP) 32

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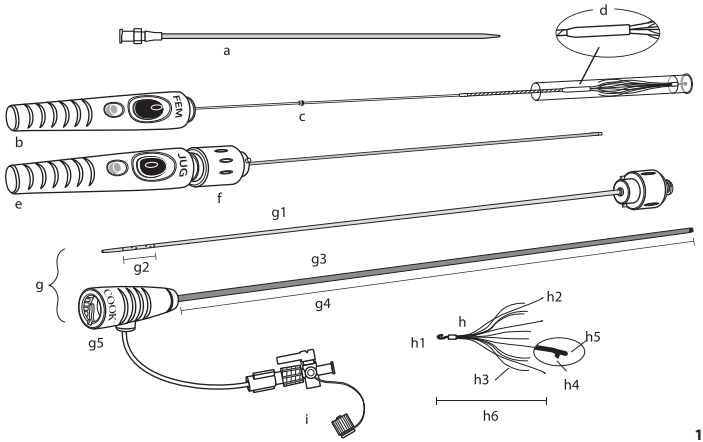
1

EN

a. Pre-dilator with hydrophilic coating
b. Femoral filter introducer with flexible tip, preloaded with the Celect Platinum filter
c. Tactile bump
d. Femoral cup
e. Jugular filter introducer with protection sheath
f. Protection sheath hub
g. Coaxial introducer system consists of:
 g1. Introducer dilator with 8 sideports and 2 radiopaque markers at the distal end
 g2. Measurement for sizing in the vena cava
 g3. Introducer sheath with radiopaque band
 g4. Sheath working length
 g5. Introducer sheath hub with Check-Flo® valve
h. Celect Platinum filter (supplied preloaded)
 h1. Hook
 h2. Primary legs
 h3. Secondary legs
 h4. Anchors
 h5. Platinum markers
 h6. Filter length
i. Three-way stopcock, plastic

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
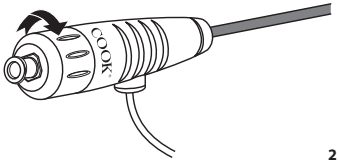
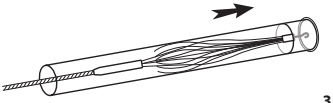
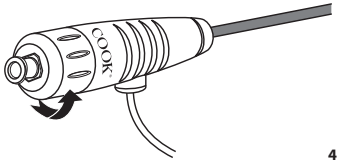
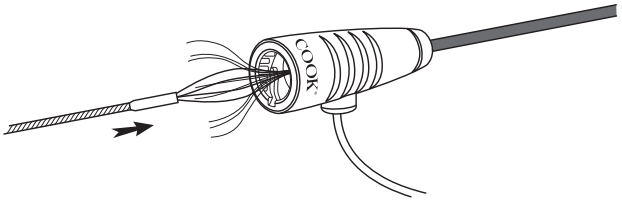
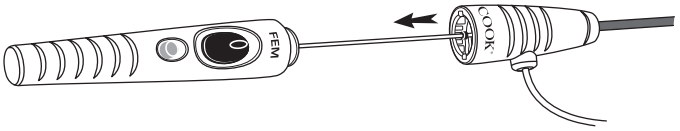
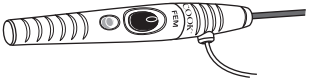
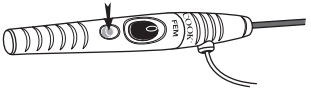
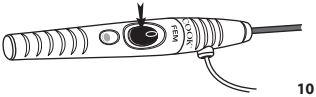
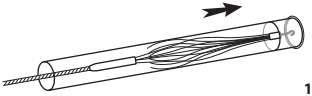
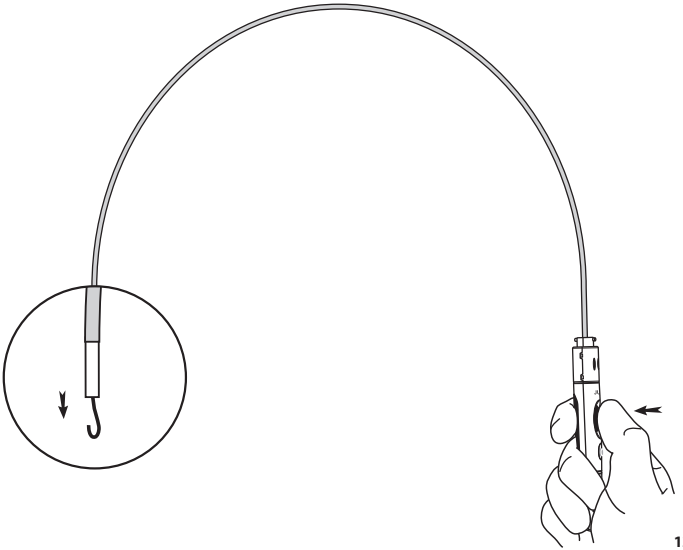
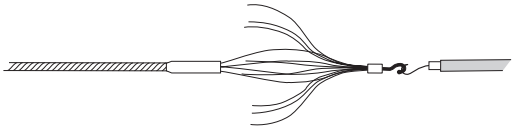
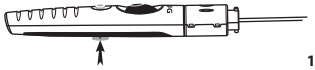
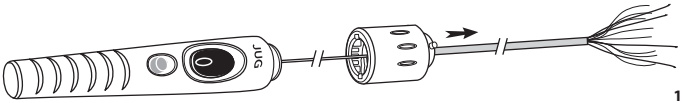
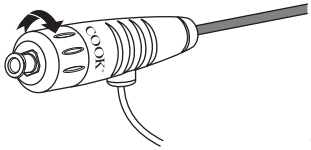
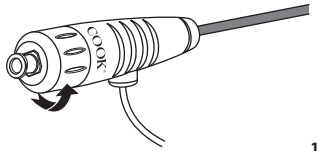
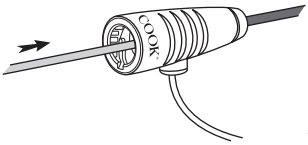
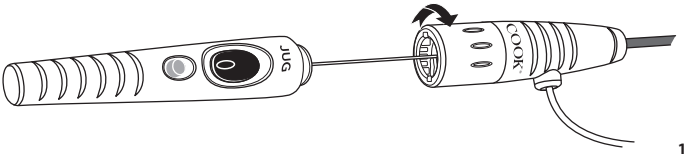
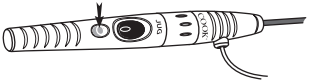
21

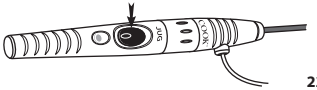22

CELECT PLATINUM® VENA CAVA FILTER SET FOR FEMORAL AND JUGULAR VEIN APPROACH

Read all instructions carefully. Failure to properly follow the information provided may lead to the device not performing as intended or injury to the patient.

STERILE—DO NOT RESTERILIZE—SINGLE USE ONLY.

1. DEVICE DESCRIPTION

1.1 Construction of Device

The Celect Platinum Vena Cava Filter Set (hereinafter referred to as the Celect Platinum filter set) consists of an inferior vena cava (IVC) filter implant and the related introducer system. The Celect Platinum Vena Cava Filter implant (hereinafter referred to as the Celect Platinum filter) is composed of a paramagnetic cobalt chromium alloy (49 mm long when compressed to a diameter of 30 mm) with platinum markers attached to the primary filter legs. The Celect Platinum filter has a total of 12 legs, four primary legs with 1 mm anchors and eight unattached secondary legs. The filter is equipped with an atraumatic, ball-tipped hook.

The Celect Platinum filter is preloaded on a femoral filter introducer with a flexible tip, which is used for placement of the filter by femoral vein approach. The Celect Platinum filter may be transferred to a jugular filter introducer also included in the universal set. The femoral and jugular filter introducers are clearly identified on their respective handles. The femoral or jugular filter introducer is used in combination with a 7 French coaxial introducer system. The coaxial introducer system includes an introducer dilator and an introducer sheath. The introducer dilator has eight sideports and two radiopaque markers 30 mm apart (end-to-end). The introducer sheath includes a radiopaque band, a hemostatic valve (Check-Flo® Valve), and a three-way stopcock in connection to the introducer sheath extension tube. When a filter introducer is used in combination with the introducer sheath the assembly will collectively be referred to as the introducer system. Finally, the Celect Platinum filter set includes a 10 French pre-dilator with hydrophilic coating for initial vessel access.

The Celect Platinum filter may be retrieved if clinically indicated; please refer to **Section 14.2.3, Optional Retrieval Procedure** for more information.

1.2 Performance Characteristics

The Celect Platinum filter is intended for percutaneous introduction into the IVC via a femoral or jugular vein access in adults. The filter is designed to capture blood clots travelling in the IVC. It is also designed to be retrieved when no longer medically indicated, although it can be permanent in some patients where retrieval may not be possible or where the benefits do not outweigh the risks of the retrieval procedure. The Celect Platinum filter is designed to ensure an appropriate radial force for a specified IVC diameter between 15 mm and 30 mm. The ball-tipped hook is used for repositioning of the filter during jugular vein delivery or retrieval of the filter. The primary and secondary legs are used to trap embolic particles. Specifically, the filter has four primary legs which are hooked outward to secure the filter to the IVC wall. The 1 mm anchors on the primary legs are designed to ensure implant fixation to the IVC. The radiopaque platinum markers at the end of the primary legs enhance implant visibility during diagnostic imaging. Likewise, the radiopaque band on the introducer sheath is intended to provide guidance during placement. The introducer sheath hemostatic valve (Check-Flo Valve) is intended to reduce potential back-flow of blood during the procedure. The stopcock is used to open, close or switch between the flow channels for system flushing. The pre-dilator is used to widen the puncture site prior to insertion of the sheath system to minimize vessel trauma. The introducer dilator is used for vena cava sizing in accordance with applicable vena cava specifications under **Section 1.6, Patient Population**, as well as injection of contrast used for visualization during diagnostic imaging.

1.3 Specifications

The filter system comes in one size. The specifications for the different parts of the system can be seen in **Table 1**.

Table 1 – Sizes of Parts Included in the Filter Set

Pre-dilator (**Fig. 1 a**)	10 Fr, 20 cm long
Measurement for sizing in the vena cava (**Fig. 1 g2**)	30 mm
Introducer sheath (**Fig. 1 g3**)	7 Fr
Sheath working length (**Fig. 1 g4**)	65 cm
Filter length when compressed to 30 mm in diameter (**Fig. 1 h6**)	49 mm

1.4 Device Compatibility

The Celect Platinum filter set is compatible with a 0.035 inch guide wire, the Cook Günther Tulip® Vena Cava Filter Retrieval Set, and the Cook CloverSnare® 4-Loop Vascular Retriever.

1.5 Qualitative and Quantitative Information

The maximum composition representative in terms of size and materials is presented in **Table 2**.

Table 2 – Qualitative and Quantitative Information

Qualitative Information Implant Material		Quantitative Information Weight (grams)
Cobalt chromium alloy*	Filter (legs and hook)	Up to 0.59
Platinum tungsten alloy	Markers on filter	Up to 0.02

*The alloy contains nickel.

1.6 Patient Population

The target population for the Celect Platinum filter is adult patients at risk of pulmonary embolism (PE) in specific situations. The Celect Platinum filter may be placed in patients requiring permanent or temporary protection to counter the risk of PE. The Celect Platinum filter is intended for percutaneous placement via a femoral or jugular vein in adults with a vena cava diameter between 15 mm and 30 mm.

The Celect Platinum filter is contraindicated for use in pregnant women and in minors/pediatric patients.

1.7 Intended User

The device is intended for use by physicians trained and experienced in diagnostic and interventional techniques.

1.8 Contact with Body Tissue

The Celect Platinum filter set is in physical contact with circulating blood in the vascular system. The contact duration for the Celect Platinum filter is long-term exposure (>30 days), while the remaining parts of the Celect Platinum filter set has a limited contact duration (≤24 hours).

1.9 Operating Principle

The Celect Platinum filter is indicated for prevention of symptomatic PE by capturing blood clots travelling in the infrarenal IVC.

The Celect Platinum filter is preloaded onto a femoral filter introducer. The filter can be transferred to the jugular filter introducer if needed. The femoral filter introducer securely fixates the primary legs and is advanced through an introducer sheath to the intended location. The jugular filter introducer securely holds the retrieval hook and is advanced through an introducer sheath to the intended location.

Fluoroscopy should be used during placement and deployment to confirm proper positioning of the Celect Platinum filter. Standard techniques for placement of vascular access sheaths, angiographic catheters and wire guides should be employed. The procedure includes use of local anesthesia, sedation and analgesia as required.

1.10 Expected Lifetime of the Device

The Celect Platinum filter is designed to be retrieved when no longer medically indicated. Current clinical practice guidelines recommend that patients with indwelling IVC filters should be routinely monitored for assessment of filter retrieval. Filter retrieval is a patient-specific, clinically complex decision that should be based on the individual patient's risk/benefit profile. Recent establishments of targeted follow-up programs, dedicated filter clinics, and institutional protocols have led to an increase in filter retrieval attempts.

In some patients, retrieval may not be possible and the device will remain implanted throughout the

life of the patient. Therefore, the expected lifetime of the IVC filters is permanent. The performance and expected lifetime of the device may be impacted by physiological processes, anatomy beyond the normal ranges, progression of an underlying vascular disease, and when the instructions for use or storage recommendations are not followed.

2. INTENDED USE

The Celect Platinum Vena Cava Filter is intended to capture blood clots traveling in the infrarenal inferior vena cava (IVC).

3. INDICATIONS FOR USE

The Celect Platinum Vena Cava Filter is indicated for prevention of symptomatic pulmonary embolism (PE) by capturing blood clots travelling in the infrarenal IVC in the following clinical situations:

- Patients with evidence of PE or deep vein thrombosis (DVT) in the IVC or lower extremity veins and either an absolute or relative contraindication to anticoagulation, a failure of anticoagulation, or high risk from complication from anticoagulation.
- Patients in the perioperative setting with a history of PE or DVT for whom anticoagulation must be interrupted.
- Patients at high risk of DVT or PE but considered unable to receive anticoagulation due to a high risk of complication.

The Celect Platinum Vena Cava Filter may be retrieved if clinically indicated; please refer to **Section 14.2.3, Optional Retrieval Procedure** for more information.

4. CLINICAL BENEFITS

For select patients who require PE prevention, IVC filtration is the only treatment option available. Therefore, PE prevention using the Celect Platinum Vena Cava Filters can be lifesaving. Analysis of the clinical data demonstrate that the incidence and type of potential adverse events related to the Celect Platinum Vena Cava Filters are acceptable and that the potential benefits to the patient population are considered to outweigh the risks. The Celect Platinum Vena Cava Filter is intended to provide the clinical benefit of prevention of PE.

5. CONTRAINDICATIONS

5.1 Filter Placement

- Megacava (diameter of the IVC > 30 mm).
- Diameter of the IVC < 15 mm.
- Extensive thrombus in the vein chosen for approach.
- Patients with risk of septic embolism.
- Use in pregnant women.
- Use in minors/pediatric patients.

5.2 Optional Filter Retrieval

- Filters with significant amounts of trapped thrombus (greater than 25% of the volume of the cone).
- Patients with ongoing high risk of PE.

6. WARNINGS

6.1 Sterile and Single Use

- Do not use the device if the sterile packaging is damaged or unintentionally opened before use. If a contaminated device is introduced to the bloodstream, it may cause an infection and lead to severe harm.
- This single-use device is not designed for reuse. If the device is reused, it could pose risks of cross contamination with microbiological agents. Attempts to reprocess (resterilize) and/or reuse may lead to device failure and/or transmission of disease.

6.2 Malfunctions/Changes in Performance

- If severe resistance is met when advancing the wire guide or the introducer system, then retract and choose a different approach. Excessive force should not be exerted.
- When power injecting contrast media, do not exceed the maximum pressure rating of 68 bar/1000 psi and flow rate of 20 mL/sec. Hand injection is also possible.
- When inserting the preloaded filter into the Check-Flo valve of the introducer sheath, hold the introducer with the flexible tip near the end, close to the filter.
- Do not attempt to rotate the preloaded filter inside the introducer system.
- Do not re-sheath the expanded filter during placement using the femoral approach.
- Do not attempt to rotate, advance, or retract the expanded filter inside the vena cava.
- Excessive force should not be exerted in placement of the filter. If deployment of the filter is not possible, it may require a replacement of the device. If a replacement of the device is not possible or if the filter does not expand correctly, it may require additional interventions or surgical removal.
- During diagnostic imaging evaluate that the filter does not show any signs of damage or defect. If the filter is damaged, it may affect the clot trapping ability of the filter or cause an obstruction of the blood flow.
- Excessive force should not be exerted to reposition or retrieve the filter as it may lead to filter breakage and/or harm to the patient. If repositioning or

retrieval of the filter is complicated, it may require additional interventions or surgical removal.

6.3 Optional Filter Retrieval

- An inferior vena cava imaging evaluation for residual captured thrombus should be performed prior to attempted retrieval.
- Never attempt to re-deploy a retrieved filter.
- For all retrievable IVC filters, retrieval becomes more challenging with time, and this is commonly due to encapsulation of the filter legs or hook (in a tilted filter) by tissue ingrowth.

7. PRECAUTIONS

7.1 General Precautions

- Standard techniques for placement of vascular access sheaths, angiographic catheters and wire guides should be employed.
- Manipulation of products (e.g., placement and retrieval) require imaging control.
- Follow the instructions thoroughly to ensure successful deployment, and to avoid any harm to the patient or damage to the device.
- Failure to store the device correctly may result in material degradation and/or damage to the device.

7.2 Malfunctions/Changes in Performance

- The Celect Platinum filter set should be used in patients with vessel diameters compatible with associated device components.
- Product (filter or introducer system) modification or alteration is not recommended, as the product's safety and effectiveness have not been established following any modifications.
- Before injecting any contrast media (by either power or hand injection) through the introducer dilator, ensure that the introducer sheath hub and introducer dilator are correctly connected.
- Filter tilt has been reported. Potential causes may include improper deployment; manipulations near an implanted filter (e.g., a surgical or endovascular procedure in the vicinity of a filter); and/or a failed retrieval attempt. Excessive filter tilt may contribute to difficult or failed retrieval, vena cava wall penetration/perforation, and/or result in loss of filter efficiency.
- Vena cava wall penetration/perforation has been reported and may be either symptomatic or asymptomatic. Potential causes may include improper deployment, and/or excessive force or manipulations near an in situ filter (e.g., a surgical or endovascular procedure in the vicinity of a filter).
- Filter fracture has been reported and may be either symptomatic or asymptomatic. Fracture of a filter leg may be due to repetitive motion on a filter leg in an unusual, stressed position, such as a filter leg penetrating/perforating the IVC or a filter leg being caught in a side branch (e.g., a renal vein). Other potential causes of filter fracture may include excessive force or manipulations near an implanted filter (e.g., a surgical or endovascular procedure in the vicinity of a filter). Retrieval of a fractured filter or filter fragments (including embolized fragments) using endovascular techniques has been reported.
- Filter or filter fragment migration and/or embolization (e.g., movement to the heart or lungs) has been reported. Filter or filter fragment movement has occurred in both the cranial and caudal direction and may be either symptomatic or asymptomatic. Potential causes may include improper deployment; deployment into thrombus; dislodgement due to large thrombus burdens; and/or excessive force or manipulations near an in situ filter (e.g., a surgical or endovascular procedure in the vicinity of a filter).
- Increased friction and/or compression at the access site during the procedure may lead to increased risk of thrombosis at the access site.
- If the introducer system or parts of the introducer system malfunctions prior to or during procedure, the device should be replaced. If the device malfunctions during the procedure, perform careful replacement to avoid injuries to the access site and vessel.

7.2.1 Malfunctions/Changes in Performance for Femoral Filter Placement

- For placement of the filter, the right femoral vein is usually preferred due to its straighter route to the IVC. The left femoral vein can be used but is more tortuous. Prior to choosing an approach, assess the patient's size and anatomy, and the location of any venous thromboses.
- The filter implant is supplied preloaded on the femoral filter introducer. Do not attempt to separate the preloaded filter introducer.
- Do not attempt to reload the filter onto the femoral filter introducer; any attempt to do so may damage the introducer and/or the filter.

7.2.2 Malfunctions/Changes in Performance for Jugular Filter Placement

- For placement of the filter, the right jugular vein is usually preferred due to its straighter route to the vena cava. An approach via the left jugular vein may be possible, depending on the patient's size and anatomy, and the location of any venous thromboses.
- The filter may be repositioned prior to final deployment by carefully advancing the introducer

sheath over the filter; repositioning the system as desired; and again, withdrawing the introducer sheath by reattaching it to the protection sheath hub, completely exposing the filter.

7.3 CMR 1A/1B and/or Endocrine Disrupting Substances

- This symbol on the label indicates that the device contains cobalt (Co) at a level above 0.1% w/w, which is a substance that is toxic to reproduction and carcinogenic (Category 1B). However, the device contains cobalt as a cobalt chromium alloy and does not release cobalt at levels which cause an increased risk of cancer or adverse reproductive effect to the patient.

7.4 Allergic Reaction

- Possible allergic reactions to cobalt, chromium, nickel, and platinum should be considered.
- Ensure that the patient does not have impaired tolerance to general, regional, or local anesthesia to avoid adverse reactions associated with the anesthetic procedure.
- Ensure that the patient is not allergic/sensitive to contrast media since the use of contrast media during the procedure and/or during postoperative imaging may cause an allergic reaction and/or other contrast-induced harms.

7.5 Training Requirements for User

The Celect Platinum filter set is intended for use by physicians trained and experienced in diagnostic and interventional endovascular techniques.

7.6 Requirements for Facilities

The device is intended for use in hospital operating rooms. Users must wear standard operating room personal protective equipment (PPE; gloves, face mask, sterile gown, etc.) and follow local guidelines for sterile procedures.

8. MRI SAFETY INFORMATION

MR Conditional

A patient with the Celect Platinum Vena Cava Filter may be safely scanned under the following conditions. Failure to follow these conditions may result in injury.

Parameter	Notes
[1] **Item Name/Identification**	Celect Platinum Vena Cava Filter
[2] **Item Manufacturer**	Cook Medical
[3] **Static Magnetic Field Strength [T]**	1.5 T or 3.0 T
[4] **Maximum Spatial Field Gradient [T/m and gauss/cm]**	20 T/m (2000 gauss/cm)
[5] **RF Excitation**	Circularly Polarized (CP)
[6] **RF Transmit Coil Type**	Whole body transmit coil, Head RF transmit-receive coil
[7] **RF Power**	Normal Operating Mode
[8] **Maximum Whole Body SAR [W/kg]**	2.0 W/kg
[9] **Scan Duration**	2.0 W/kg whole body average SAR for 15 minutes of continuous RF (a sequence or back to back series/scan without breaks). Under the scan conditions defined above, the Celect Platinum Vena Cava Filter is expected to produce a maximum temperature rise of less than 4.2 °C after 15 minutes of continuous scanning.
[10] **MR Image Artifact**	The presence of this implant may produce an image artifact of 24 mm.
If information about a specific parameter is not included, there are no conditions associated with that parameter.	
Follow the MRI safety information to avoid excessive heating, torque, and/or deflection, which may cause injury to the vessel.	
Image artifacts may occur, which may prolong diagnostic time and/or require additional imaging.	

9. POTENTIAL ADVERSE EVENTS

Potential adverse events that may occur include the following:

- Access site thrombosis/occlusion
- Air embolism
- Arrythmia
- Back or abdominal pain
- Blood loss
- Branch vessel occlusion
- Cardiac damage
- Cardiac tamponade
- Damage to the vena cava
- Death
- Deep vein thrombosis
- Edema
- Extravasation of contrast material
- Failure of filter expansion/incomplete expansion
- Filter fracture
- Filter malposition
- Filter migration
- Filter or filter fragment embolization
- Hematoma at vascular access site
- Hemorrhage
- Infection at vascular access site
- Intimal tear
- Obstruction of blood flow
- Pneumothorax
- Postphlebitic syndrome
- Pulmonary embolism
- Retrieval failure
- Trauma to adjacent structures
- Unacceptable filter tilt
- Vascular trauma
- Vena cava occlusion or thrombosis
- Vena cava penetration
- Vena cava perforation
- Vena cava stenosis

10. STERILIZATION

All devices are sterilized using ethylene oxide (EO) gas.

11. HOW SUPPLIED

This device is accompanied by an implant card, that should be given to the patient after it has been completed by the healthcare professional.

Final devices are supplied in peel-open packages. Keep the device dry and away from sunlight. Do not use after the expiration date printed on the label.

12. INSPECTION OF DEVICE

Inspect the device thoroughly including all levels of the packaging to verify that there is no damage prior to use. Inspect and confirm that the sterile barrier has not been compromised in any way. Inspect and confirm that the device corresponds to the label and the IFU.

Do not use the device if the sterile packaging is damaged or unintentionally opened before use.

13. DEVICE PREPARATION

13.1 Device Selection

The Celect Platinum filter comes in one size and is available with a femoral or jugular introducer system.

14. INSTRUCTIONS FOR USE

14.1 Patient Preparation

The physician should follow institutional protocols and local guidelines.

14.2 Step-by-Step Guide

The step-by-step guide is divided in descriptions of femoral and jugular approach.

14.2.1 Femoral Approach

14.2.1.1 Preparation

1. Flush the introducer sheath and the introducer dilator.
2. Advance the introducer dilator through the middle of the Check-Flo valve on the introducer sheath. Secure the introducer dilator to the introducer sheath by twisting the dilator hub clockwise until a click is felt. (**Fig. 2**)
3. Remove the filter protection tube. (**Fig. 3**)

14.2.1.2 Filter Placement

4. Access the chosen femoral vein using the Seldinger technique.
5. Perform diagnostic imaging to confirm a single IVC, measure the IVC diameter, check for thrombus, and establish the position of the renal veins.
6. Place a supportive 0.035 inch wire guide in the IVC.
7. If necessary, dilate the puncture site with the 10 French pre-dilator.
8. Remove the pre-dilator and advance the coaxial introducer system over the wire guide until the tip of the introducer sheath lies approximately 1 cm caudal to the lowest renal vein.
9. Remove the wire guide.
10. Perform diagnostic imaging to verify the position of the introducer sheath tip (or radiopaque marker) approximately 1 cm caudal to the lowest renal vein.

CAUTION: Before injecting contrast media by either power or hand injection through the introducer dilator, ensure that the introducer sheath hub and introducer dilator are correctly connected.

WARNING: When using a power injector, do not exceed the maximum pressure rating of 68 bar/1000 psi and flow rate of 20 mL/sec.

11. When correct position is established, twist the introducer dilator hub counterclockwise and remove the introducer dilator. (**Fig. 4**)
12. Place the preloaded filter into the Check-Flo valve of the introducer sheath. (**Fig. 5**)

 WARNING: Hold the femoral filter introducer near the end, close to the filter, to avoid kinking the flexible tip.
13. Advance the femoral filter introducer until the Check-Flo valve contacts the tactile bump on the femoral filter introducer. This will place the hook of the filter inside the introducer sheath at the radiopaque band. Verify that the position of the hook is inside the introducer sheath and still caudal to the renal veins.

 WARNING: Do not rotate the preloaded filter inside the introducer system.

 WARNING: Do not exert excessive force to advance the filter through the introducer system.
14. Stabilize the introducer system, withdraw the introducer sheath, (**Fig. 6**) and connect it to the handle of the femoral filter introducer. (**Fig. 7**) At this point the filter is expanded, still connected to the femoral filter introducer. (**Fig. 8**)

 CAUTION: Attempting to retract the filter at this point of the deployment sequence could damage the secondary legs or caval wall.
15. Proper position can now be verified by diagnostic imaging.

 WARNING: Do not rotate the expanded filter inside the vena cava. Doing so may compromise the performance of the filter.

 CAUTION: Injection of contrast media must not be performed unless the femoral cup (indicated as position d in Fig. 1) is completely free of the introducer sheath. Use the radiopaque band for positioning.
16. Verify that the introducer sheath hub and femoral filter introducer handle are connected to ensure that the femoral cup is completely free of the introducer sheath before filter release.
17. When the filter position is correct, push the red safety button to prepare filter release. (**Fig. 9**)
18. Push the release button completely to ensure proper release of the filter. (**Fig. 10**) Repositioning of the filter is no longer possible. The filter is now released.
19. Perform diagnostic imaging to verify filter position.

 NOTE: Institutional standard of care should be followed for removing the introducer sheath and providing hemostasis to prevent bleeding at the vascular access site.

14.2.2 Jugular Approach

14.2.2.1 Preparation

To prepare for jugular approach, the filter must be transferred from the femoral filter introducer to the jugular filter introducer.

1. Remove the protection tube from the filter on the femoral filter introducer. (**Fig. 11**)
2. Push and hold the release button on the jugular filter introducer handle to advance the grasping hook beyond the protection sheath. While holding the jugular filter introducer in a soft bend according to illustration, (**Fig. 12**) catch the hook of the preloaded filter on the femoral filter introducer and release the button to firmly grasp the filter. (**Fig. 13**)

 NOTE: Do not kink the jugular filter introducer while bending it.
3. Secure the lock on the jugular filter introducer handle by pushing the red knob on the back side. (**Fig. 14**)
4. To release the filter from the femoral filter introducer, push the red safety button, followed by the release button on the femoral filter introducer. (**Fig. 9** and **10**)
5. Straighten the system and then advance the protection sheath hub on the jugular filter introducer until a confirmed stop is felt (approximately 6 cm). (**Fig. 15**) This ensures that the filter is completely inside the tip of the protection sheath.

 The product is now prepared for use.
6. Flush the introducer sheath and the introducer dilator.
7. Advance the introducer dilator through the middle of the Check-Flo valve on the introducer sheath. Secure the introducer dilator to the introducer sheath by twisting the dilator hub clockwise until a click is felt. (**Fig. 16**)

14.2.2.2 Filter Placement

8. Access the chosen jugular vein using the Seldinger technique.
9. Perform diagnostic imaging to confirm a single IVC, measure the IVC diameter, check for thrombus, and establish the position of the renal veins.

10. Place a supportive 0.035 inch wire guide in the IVC.
11. If necessary, dilate the puncture site with the 10 French pre-dilator.
12. Remove the pre-dilator and advance the coaxial introducer system over the wire guide until the tip of the introducer sheath lies approximately 5 cm caudal to the lowest renal vein.
13. Remove the wire guide.
14. Perform diagnostic imaging to verify position of the introducer sheath tip (or radiopaque marker) approximately 5 cm caudal to the lowest renal vein.

 CAUTION: Before injecting contrast media by either power or hand injection through the introducer dilator, ensure that the introducer sheath hub and introducer dilator are correctly connected.

 WARNING: When using a power injector, do not exceed the maximum pressure rating of 68 bar/1000 psi and flow rate of 20 mL/sec.
15. When correct position is established, twist the introducer dilator hub counterclockwise and remove the introducer dilator. (**Fig. 17**)
16. Place the jugular filter introducer, with the protection sheath containing the preloaded filter, into the Check-Flo valve of the introducer sheath. Advance the jugular filter introducer with the protection sheath into the introducer sheath. (**Fig. 18**)

 WARNING: Do not rotate the preloaded filter inside the introducer system.
17. Connect the introducer sheath hub and protection sheath by twisting clockwise until a click is felt. (**Fig. 19**) The filter is now positioned at the radiopaque band of the introducer sheath. The hook of the filter should be caudal to the renal veins.

 WARNING: Do not exert excessive force to advance the filter through the introducer sheath.
18. Stabilize the introducer system and withdraw the introducer sheath and protection sheath until the protection sheath and jugular filter introducer handle are in contact with one another. At this point the filter is expanded, still connected to the jugular filter introducer. (**Fig. 20**)

 WARNING: Do not rotate the expanded filter inside the vena cava. Doing so may compromise the performance of the filter.
19. If the filter is not in the desired position, carefully advance the introducer sheath over the filter. Reposition the system as desired, and again withdraw the introducer sheath and protection sheath until the protection sheath and jugular filter introducer handle are in contact with one another, completely exposing the filter.
20. When the filter position is correct, push the red safety button to prepare filter release. (**Fig. 21**)
21. Push the release button completely to ensure proper release of the filter. (**Fig. 22**) Repositioning of the filter is no longer possible. The filter is now released.

 NOTE: Excessive tension during deployment may prevent the filter from releasing when the release mechanism is activated.
22. Perform diagnostic imaging to verify filter position.

 NOTE: Institutional standard of care should be followed for removing the introducer sheath and providing hemostasis to prevent bleeding at the vascular access site.

14.2.3 Optional Retrieval Procedure

The Celect Platinum filter implant may be retrieved. The filter was designed to be retrieved with the Günther Tulip Vena Cava Filter Retrieval Set. It may also be retrieved with the CloverSnare 4-Loop Vascular Retriever. Please refer to the Instructions for Use provided with the Günther Tulip Vena Cava Filter Retrieval Set or the CloverSnare 4-Loop Vascular Retriever (not included in the filter set) for additional information on the optional retrieval procedure.

15. REQUIRED EQUIPMENT AND ACCESSORIES

15.1 Materials Required

- Wire guide
- Needle
- Catheter
- Introducer
- Fluoroscope with digital angiography capabilities (C-arm or fixed unit)
- Contrast media
- Power injector
- Syringe
- Flushing solution
- Sterile gauze pads

15.2 Materials Recommended

- Micropuncture access set
- In case of misplacement or inadvertent detachment of the Celect Platinum filter the Günther Tulip® Vena Cava Filter Retrieval Set or the CloverSnare® 4-Loop Vascular Retriever can be used for filter retrieval.

16. DISPOSAL OF DEVICE

After the procedure, this device may be contaminated with potentially infectious substances of human

origin and should be disposed of in accordance with institutional guidelines.

17. REFERENCES

These instructions for use are based on experience from physicians and/or their published literature. Refer to your local Cook Medical sales representative for information on available literature.

18. PATIENT COUNSELING

Please inform the patient as necessary of the relevant warnings, precautions, contraindications, measures to be taken and limitations of use that the patient should be aware of.

Furthermore, the patient should be informed about the residual risks when treated with this device and the potential adverse events related to the use of this device.

The physician must complete the implant card(s) for the patient. The patient should always carry it with them and refer to the card any time they visit additional health care practitioners, particularly for any additional diagnostic procedures (e.g., MRI).

18.1 Symptoms

The physician must advise the patient to contact their healthcare provider if he/she experiences changes in their medical condition or in case of concern.

18.2 Follow Up Information for the Patient

Follow-up is patient-specific and clinically complex; therefore, the physician must make an individual follow-up schedule for the patient.

18.3 Patient Implant Information

This information will be available for the patient in layman language (except for the MRI information, which is aimed at healthcare professionals) through cookmedical.eu/patient-implant-information. The physician must ensure that the patient has received the implant card giving access to the digital patient facing information. If the patient is not able to access the electronic information, the physician must provide the patient with a printed copy.

18.4 Summary of Safety and Clinical Performance (SSCP)

You can use the basic unique device identifier (BUDI) for this device (0827002-WCE-MDF0000624-GN) to find the Summary of Safety and Clinical Performance (SSCP) on the EUDAMED website. When the EUDAMED is available, use the following link: ec.europa.eu/tools/eudamed.

19. SERIOUS INCIDENT REPORTING

If any serious incident has occurred in relation to the device this should be reported to Cook Medical, and this shall also be reported to the competent authority the device was used in.

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© COOK 2024
William A. Cook Australia Pty. Ltd.